




Complex Rehab Systems

Quick Start Guide

Version 20.4

1. Log Into the CRS App	2
2. Client Records	3
3. Evaluations	5
4. Scanning	6
A. Scan a Person or Object	6
B. Take 3D Snapshots	11
5. Measurements	14
6. Summary Report	17
7. Technical Support	18

1. Log Into the CRS App

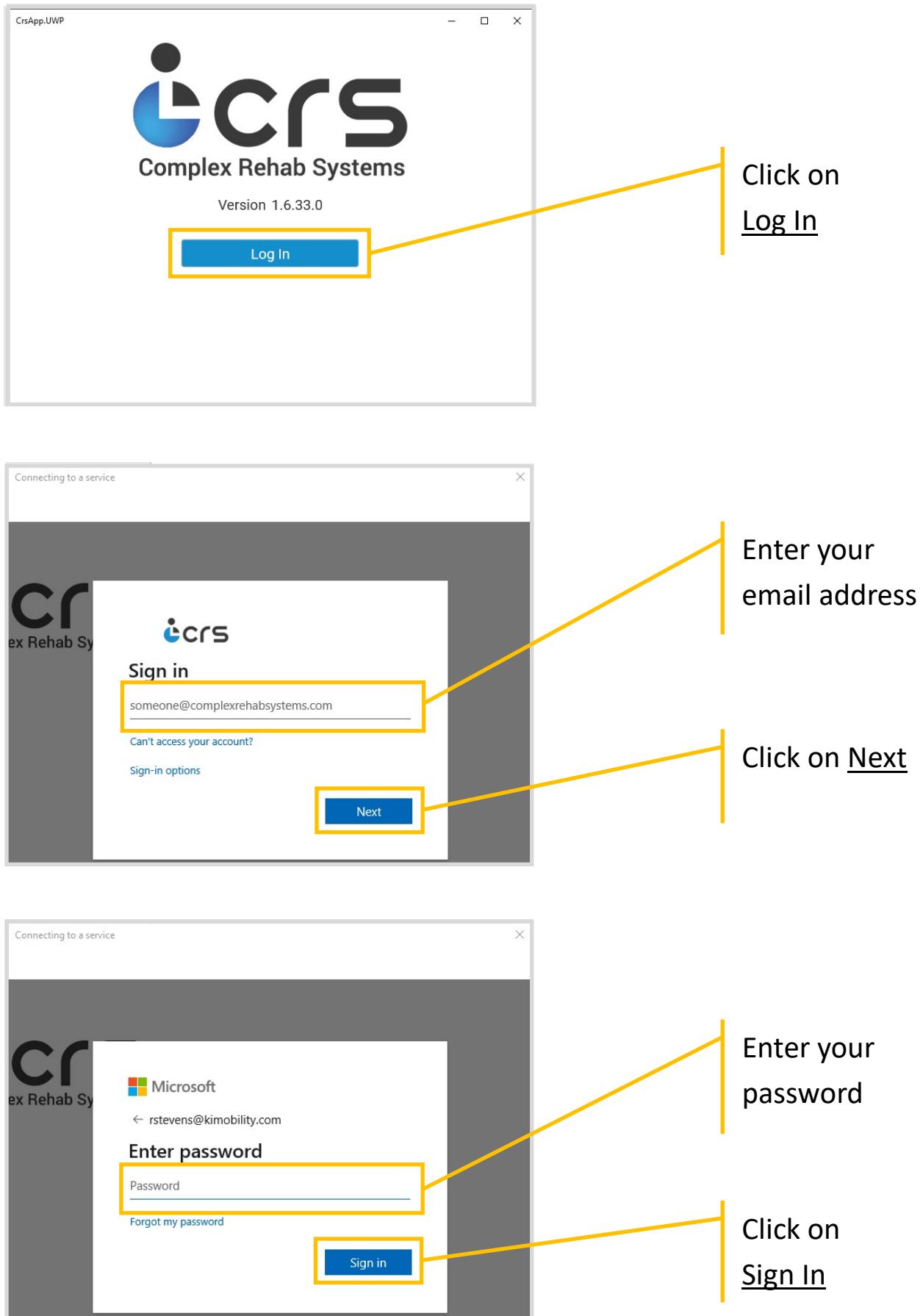

2. Client Records

Clients

Order by: Age Branch: Enter a Branch Clients: All clients

NAME	UPDATED	BRANCH	EXT. ID	STATUS	SYNC
Jones, John	03/10/2020	1		Unsynchronized Changes	
NAME ↑	UPDATED ↑	BRANCH	DOB	AGE ↑	MY CLIENTS
Stevens, Renee Elizabeth	03/06/2020		01/01/2000	20	
Saulsberry, Mitch	02/28/2020		11/15/1996	23	
Jones, Jane	03/10/2020	1	01/06/1986	34	
Jones, John	03/10/2020	1	01/06/1985	35	
Smith, Jane	12/13/2016	1	01/06/1983	37	
Smith, John	10/11/2016	1	01/06/1982	38	
Budd, Christopher	03/11/2020		12/16/1981	38	
Doe, Jane	06/07/2016	1	01/06/1979	41	
Smith, Joe	01/17/2020		01/17/1973	47	

Add a new Client

Or select an existing Client

Clients / Jones, John

Created on: 03/16/2016 | Updated: 03/10/2020

01/06/1985 | 35 years old | SSN: | Branch: 1

EVALUATION	CREATED	VISITED
03/10/2020-14:55:55	03/06/2020	03/06/2020
03/05/2020-12:06:29	03/04/2020	03/04/2020

Create a New Evaluation

Select an existing Evaluation

2. Client Records, continued . . .

Add Client

First Name		
Middle Name		
Last Name		
SSN	XXX-XX-XXXX (OPTIONAL)	
Gender	<input type="button" value="Male"/>	<input type="button" value="Female"/>
Date of Birth	January	1
Branch	Enter a Branch	
<input type="button" value="Cancel"/> <input style="outline: 2px solid yellow; border: 1px solid yellow;" type="button" value="Review"/>		

For New Clients Only

(For Existing Clients, proceed to the next page)

Enter Client information

And click on Review

Review Details

CLIENT SELECTED	
Create Client	
FULL NAME	Jason Z Kline
GENDER	Male
DATE OF BIRTH	02/01/1970

Review client information

IF ALL IS CORRECT,
Click on Confirm

3. Evaluations

Clients / Jones, John

01/06/1985 | 35 years old | SSN: | Branch: 1

EVALUATION	CREATED	VISITED
03/10/2020-14:55:55	03/06/2020	03/06/2020
03/05/2020-12:06:29	03/04/2020	03/04/2020

Created on: 03/16/2016 | Updated: 03/10/2020

New Evaluation

Click on an Existing Evaluation
Or click on New Evaluation

Evaluation Homepage

Clients / Jones, John / Evaluation

01/06/1985 | 35 years old | SSN: | Branch: 1

Work ID: Visited At: 03/11/2020 16:54:00

- Clinical Evaluations
- Notes
- Scanning
- Measurements
- Attachments
- Scratch Pad
- Summary

Contact Information

Physicians/Therapists

Referral

Goals

Functional & Sensory Abilities

ADL Status

Add Collaborator

Generate Report

Refresh

Save

Evaluation Sections
Clinical Information Cards
Scroll down for more Client Information

4. Scanning

Clients / Jones, John / Evaluation

Created on 03/11/2020 | Updated on 03/11/2020

01/06/1985 | 35 years old | SSN: | Branch: 1

Work ID: [] Visited At: 03/11/2020 16:54:00 []

Scanning

Contact Information

Physicians/Terapists

Referral

Goals

Functional & Sensory Abilities

ADL Status

Add Collaborator

Generate Report

Refresh

Save

Favorites All

Click on the
Scanning
Feature

Clients / Jones, Jane / Evaluation

Created on 03/05/2020 | Updated on 03/10/2020

01/06/1986 | 34 years old | SSN: | Branch: 1

Work ID: [] Visited At: 03/05/2020 08:53:00 []

SCAN	UPDATED	CREATED
Scan 8	03/05/2020	03/05/2020
Scan 7	03/05/2020	03/05/2020
Scan 6	03/05/2020	03/05/2020
Scan 5	03/05/2020	03/05/2020
Scan 4	03/05/2020	03/05/2020
Scan 3	03/05/2020	03/05/2020
Scan 2	03/05/2020	03/05/2020
Scan 1	03/05/2020	03/05/2020

New Scan

Import Scan


Refresh

Save

Create a
New Scan


View an
Existing Scan

These are your scanning choices:

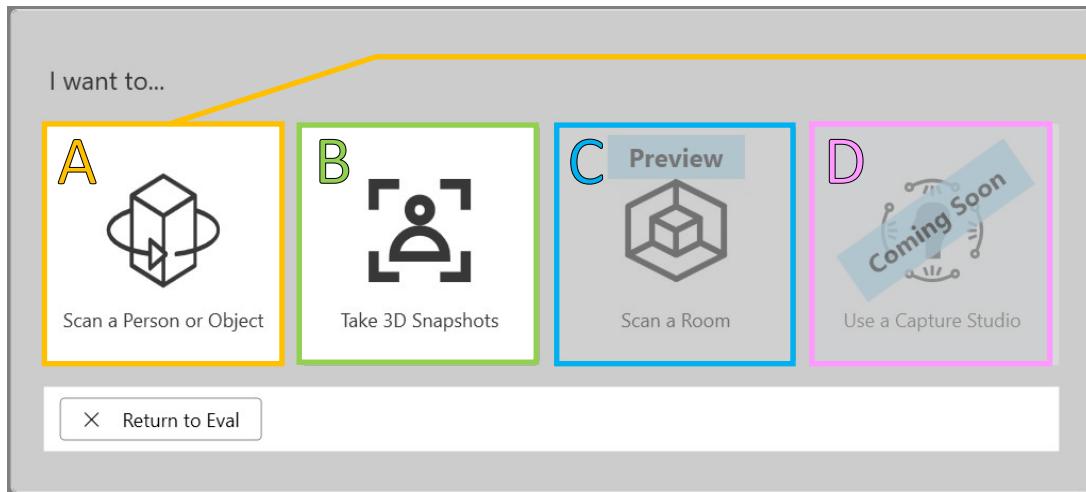

Click on Scan a Person or Object

A. Scan a Person or Object

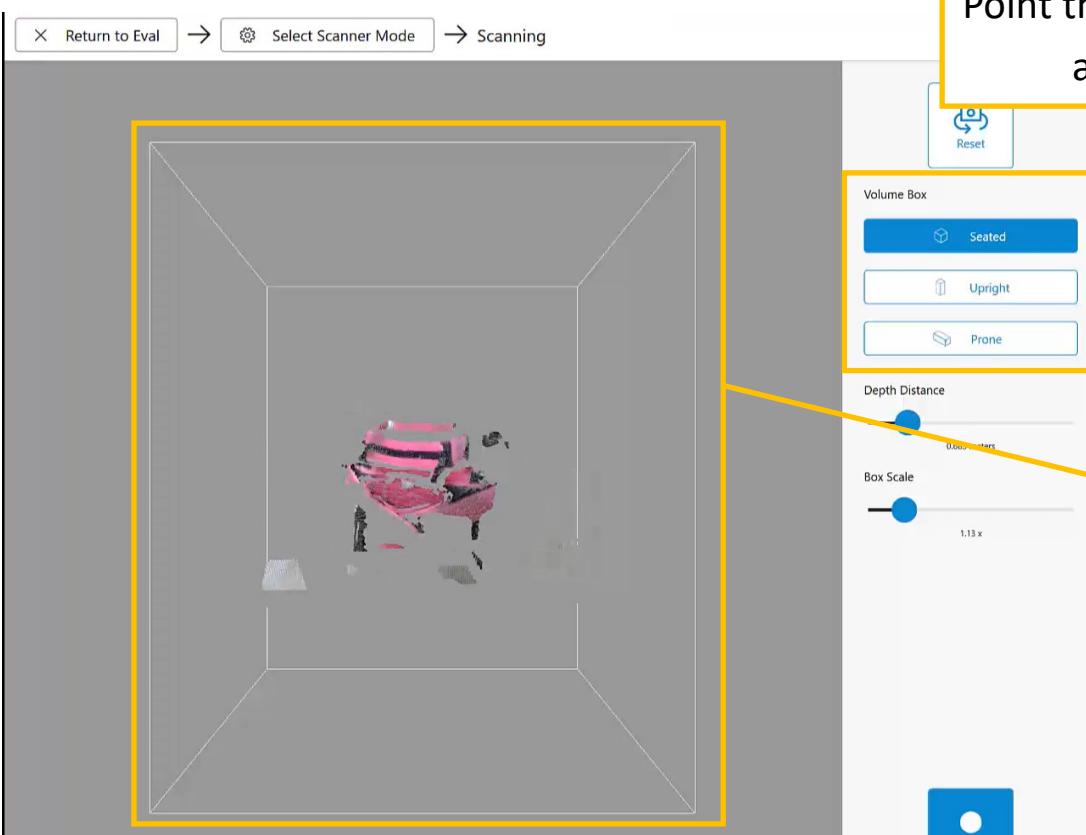


Point the scanning camera at the subject

Select the position of the subject

Touch the screen, click and drag to rotate the 3D Bounding Box.


Move the camera closer or further from the subject until they are seen in the center of the 3D Bounding Box.

4A. Scan a Person or Object, continued . . .

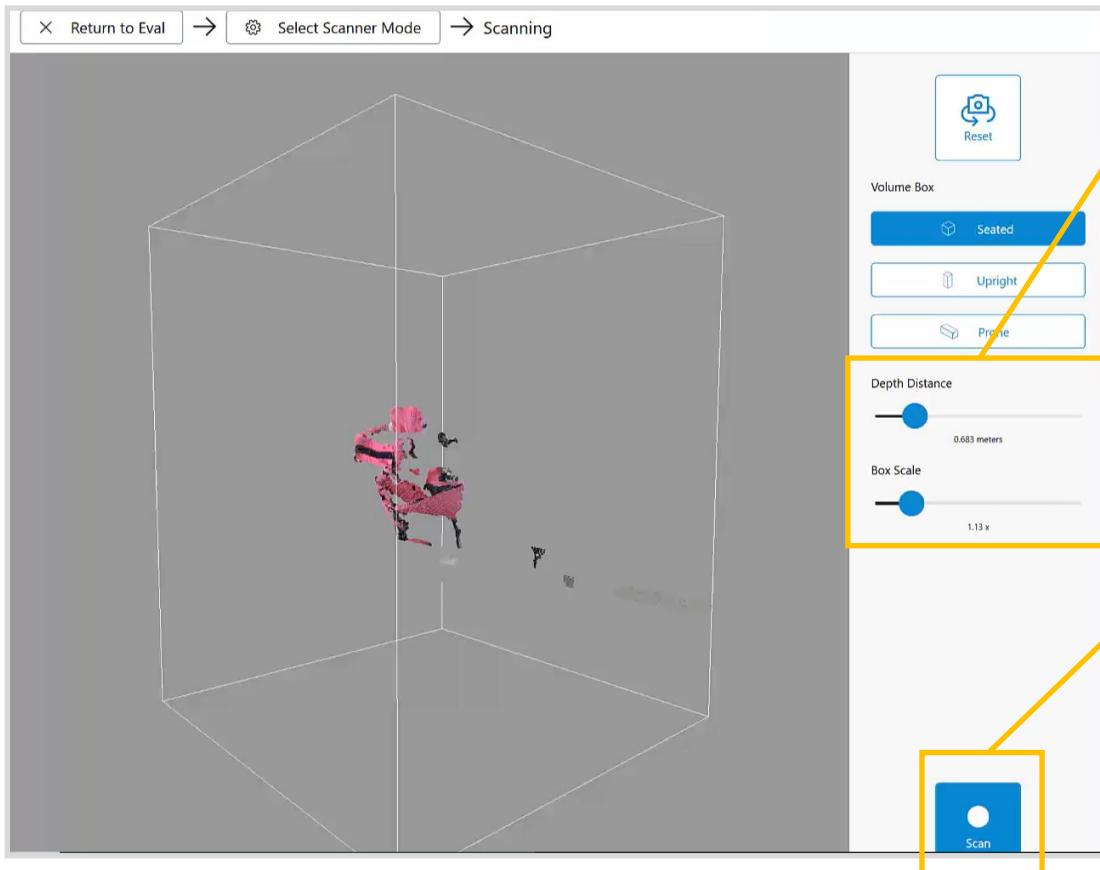

You can adjust the Depth Distance and Box Scale to help center the subject

Click on Scan to begin

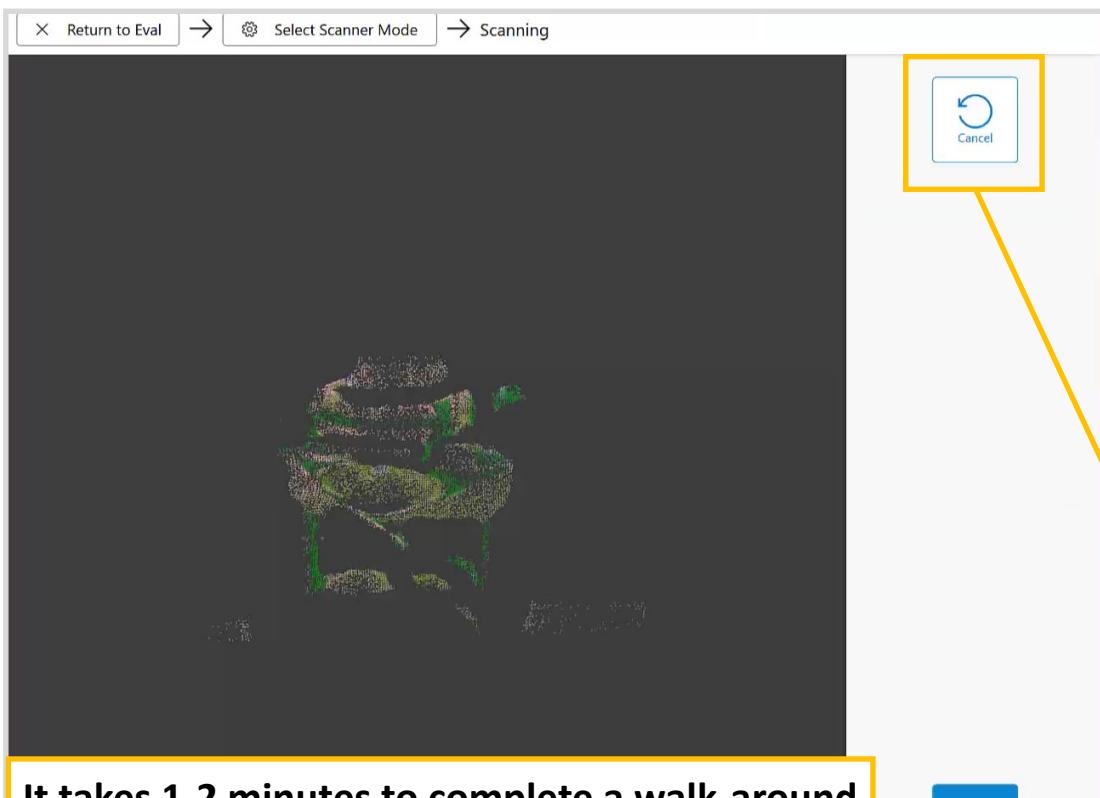


Maintaining the same distance, slowly walk around the subject, moving the scanner up and down.


It takes 1-2 minutes to complete a walk-around scan of a seated subject.

If you encounter an error or problem, you can click on Cancel and restart the scan

4A. Scan a Person or Object continued . . .


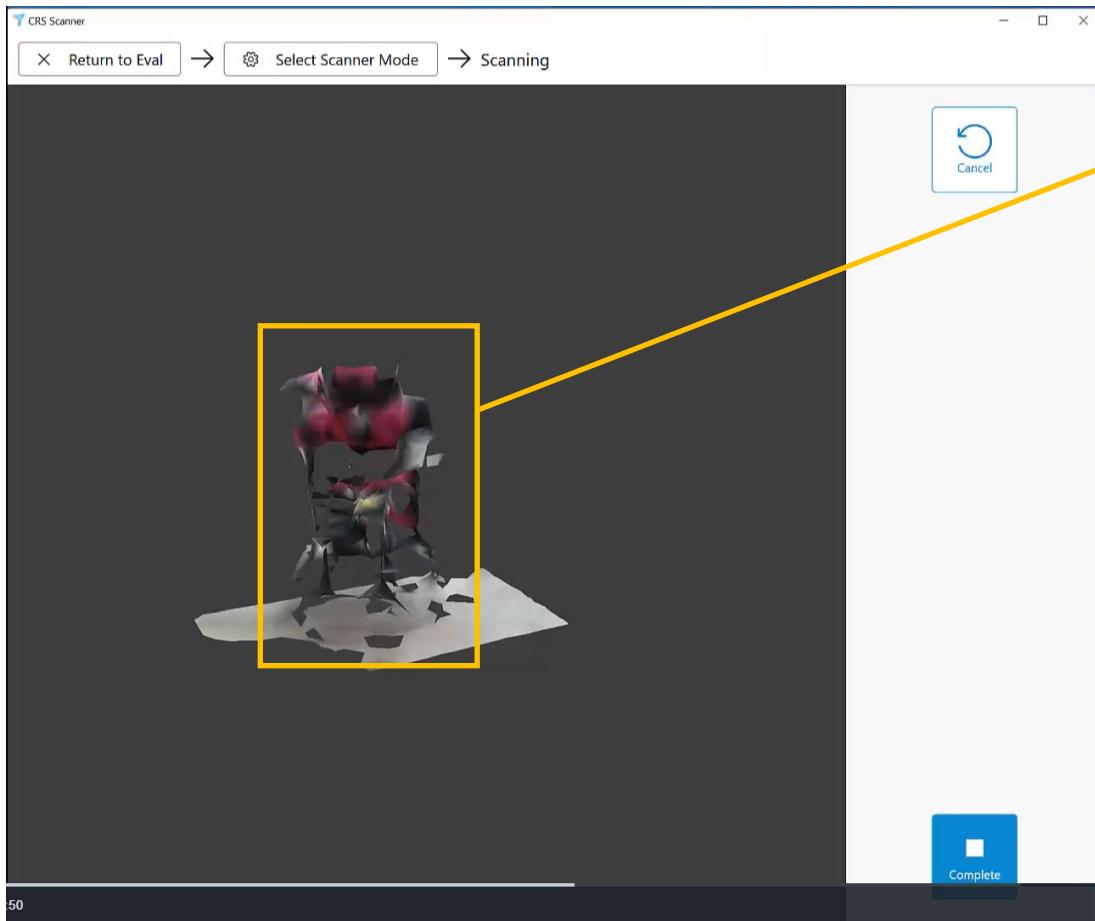


Red shadows indicate a Scanner Misalignment. Try to match the shadow with the scanned location to return to regular scanning.

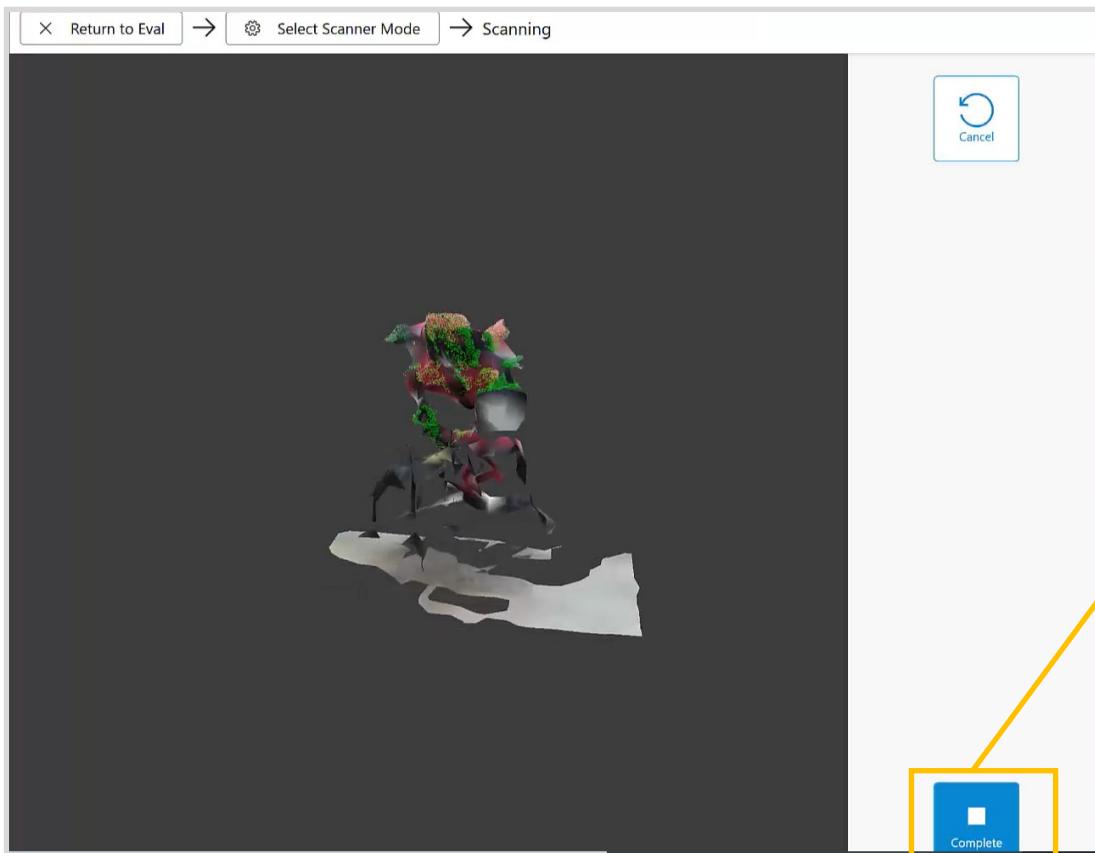


When done,
click
Complete

Review your Scan

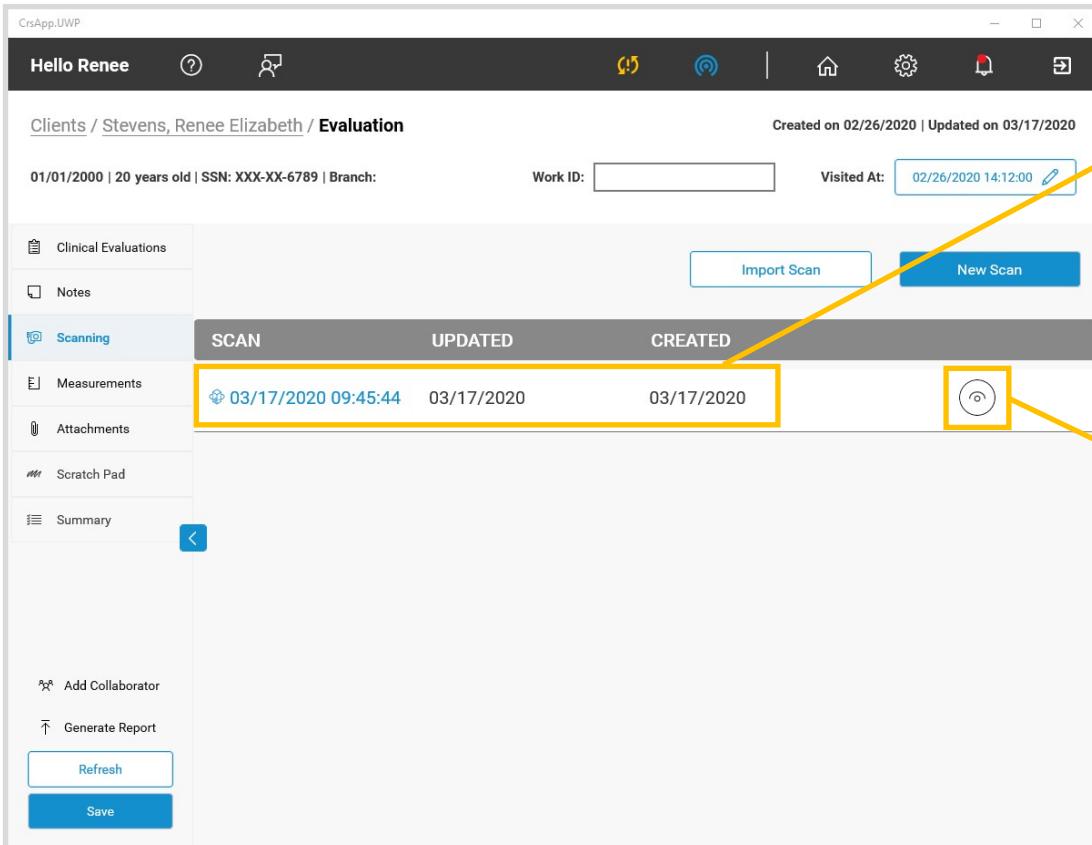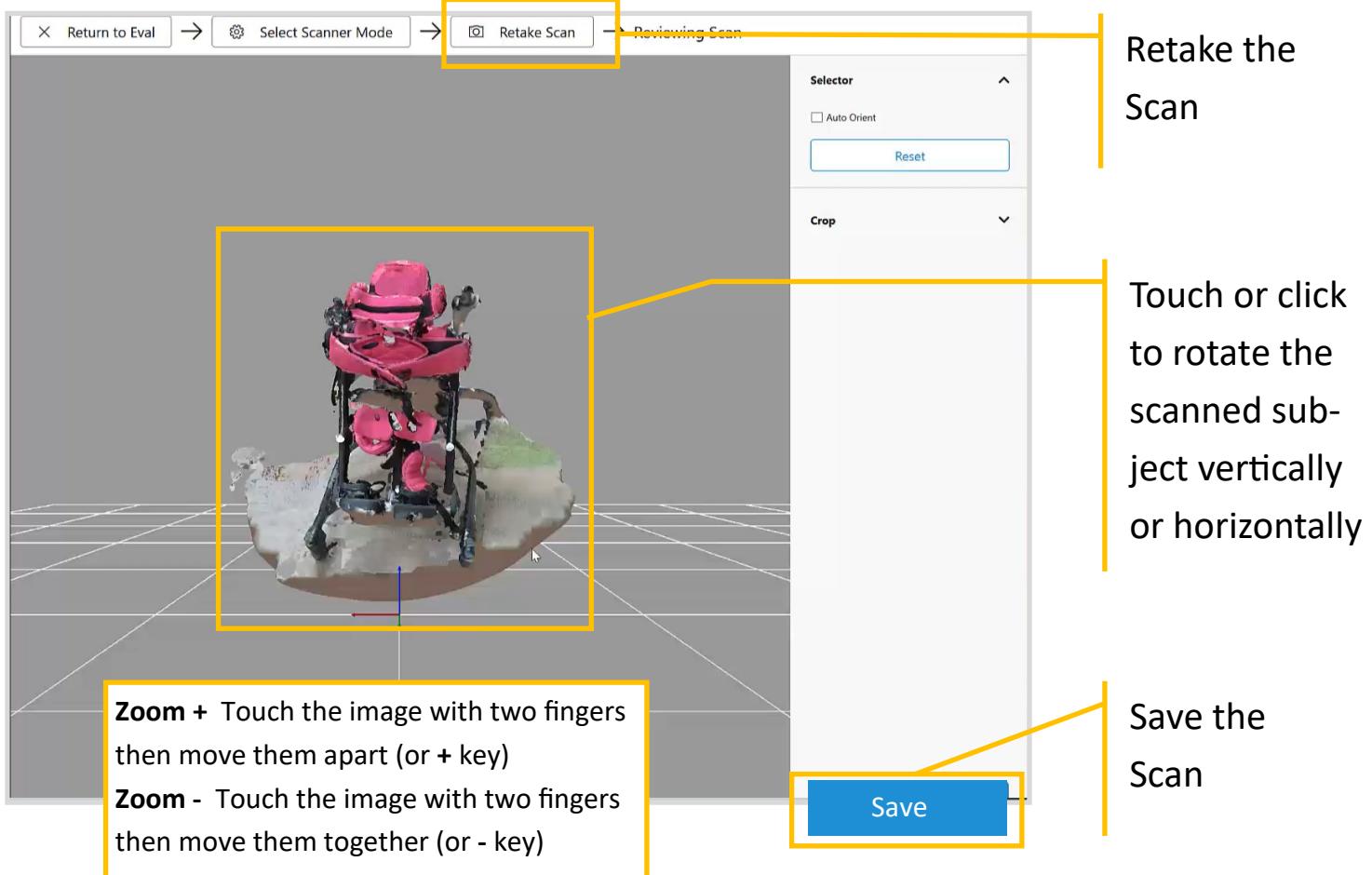


4B. Take 3D Snapshots



From the Evaluation, click on the Scanning tab and then on New Scan.

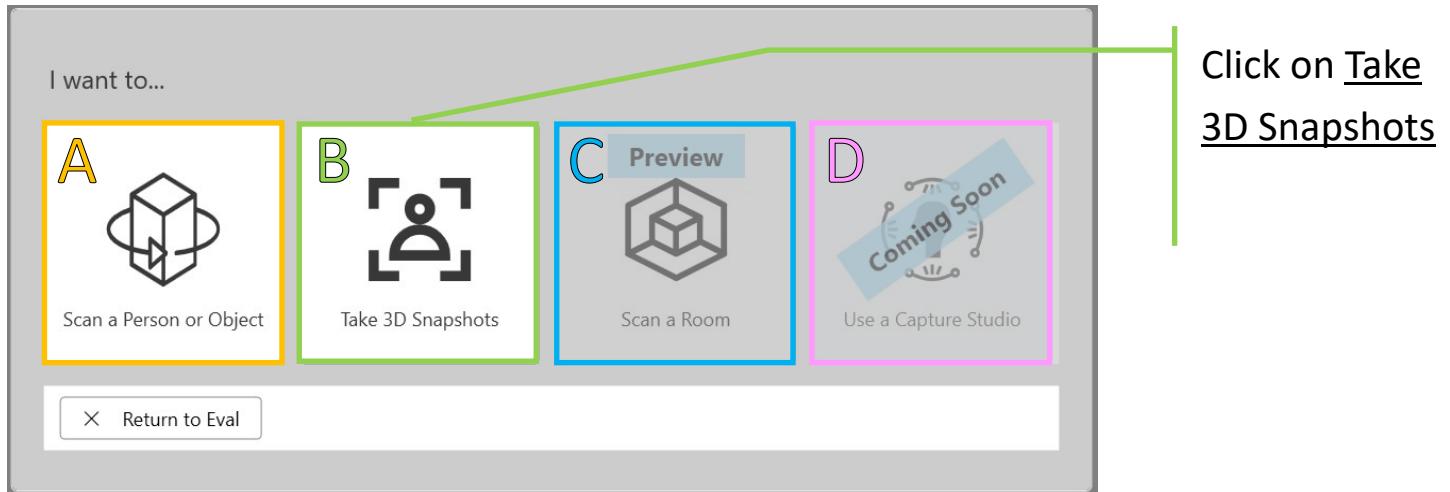


Point your camera at the client.

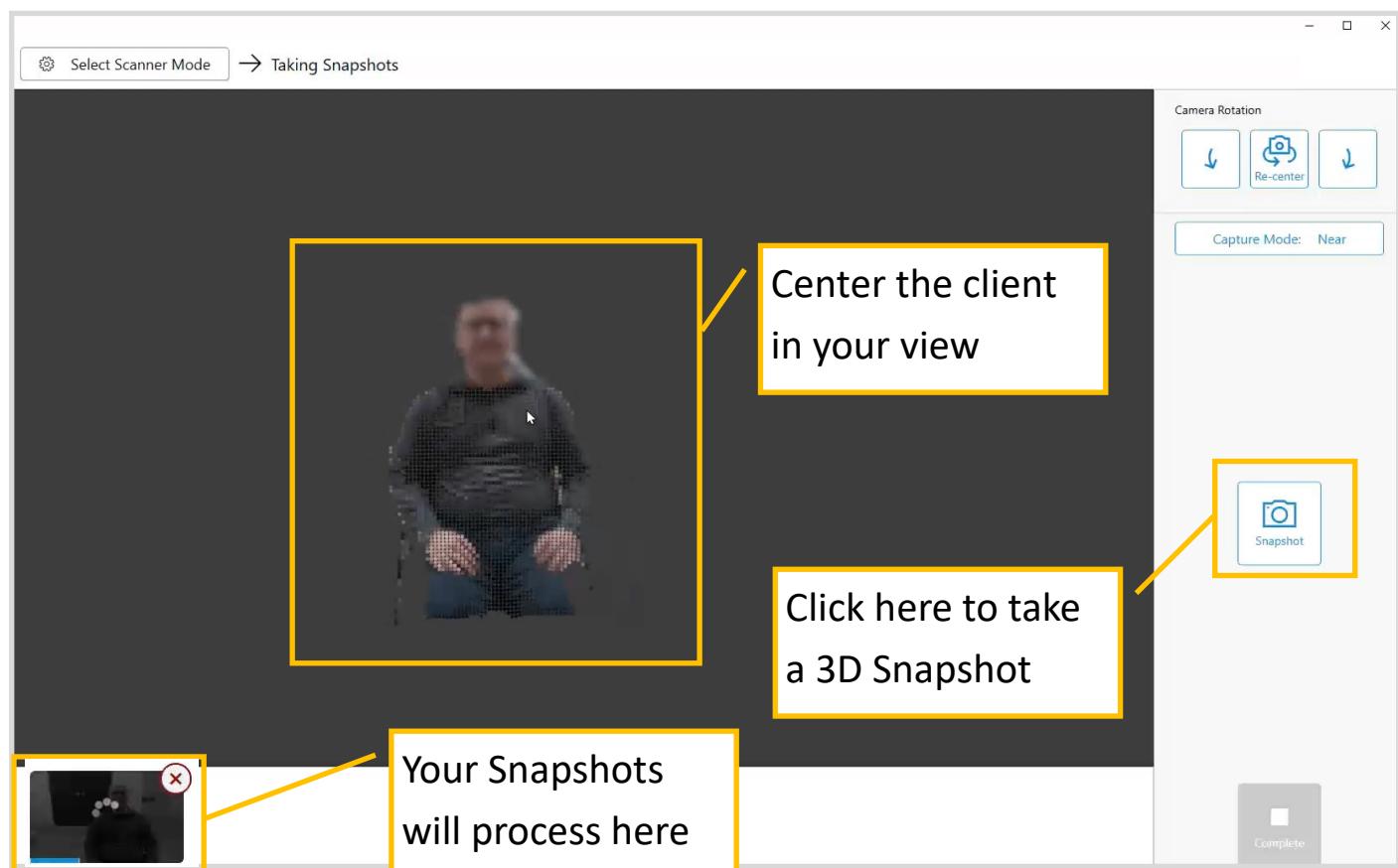

4B. Take 3D Snapshots, continued . . .

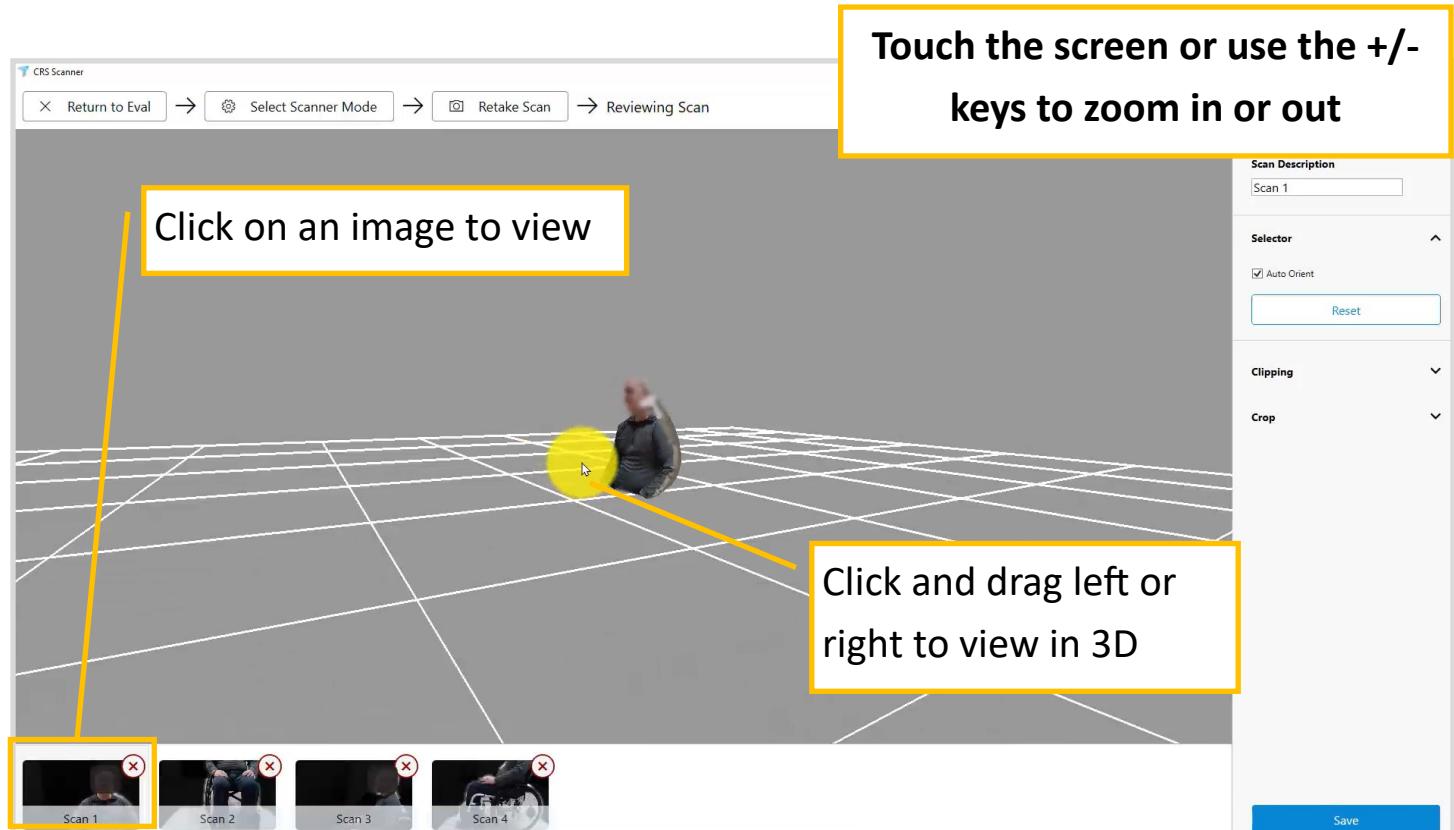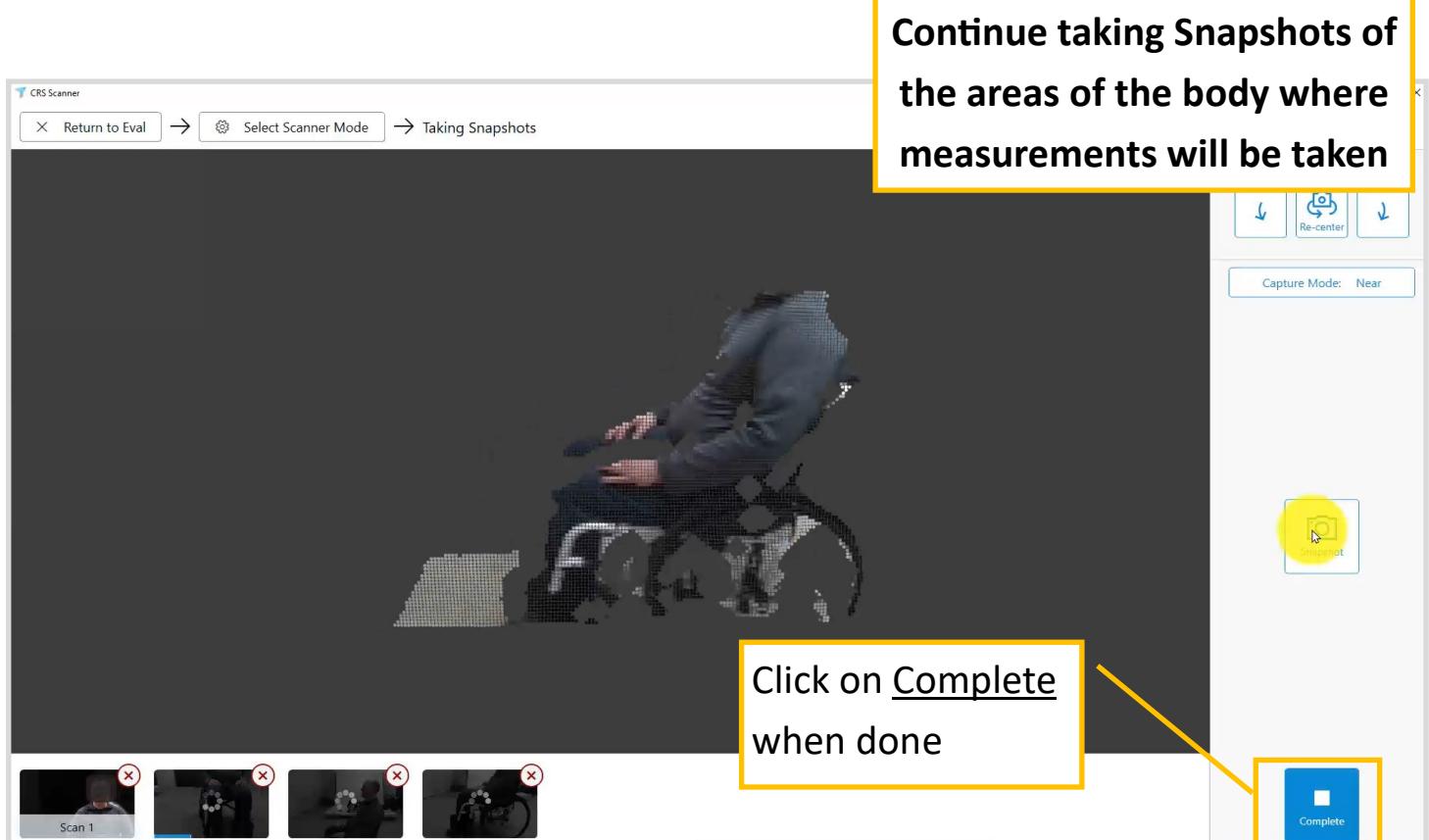

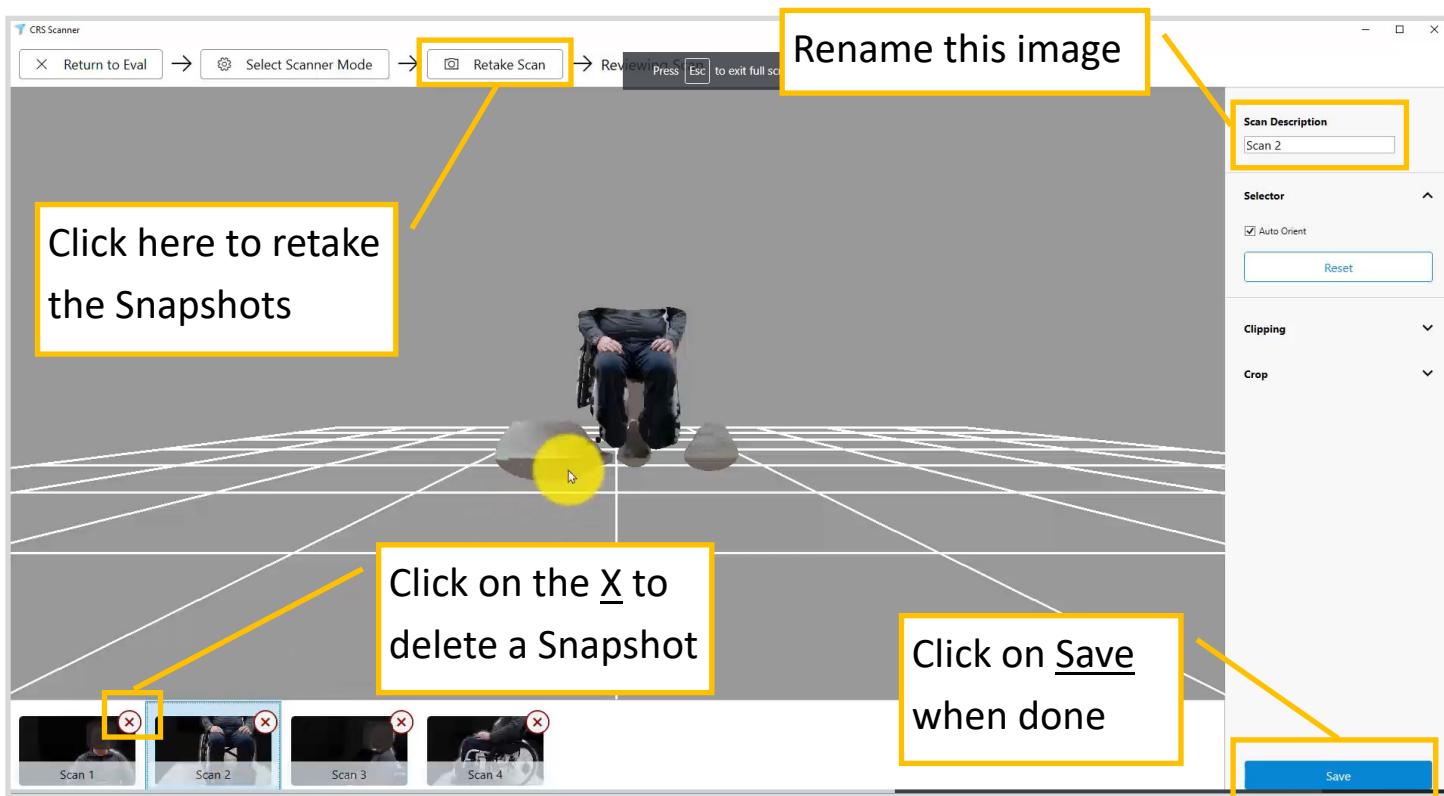

The screenshot shows the CrsApp.UWP application. On the left, there's a sidebar with navigation links: 'Clinical Evaluations', 'Notes', 'Scanning' (which is currently selected), 'Measurements', 'Attachments', 'Scratch Pad', and 'Summary'. The main area displays a table of snapshots. The table has columns for 'SCAN', 'UPDATED', and 'CREATED'. Three rows are visible, each representing a scan: 'scan 2' (updated and created on 04/08/2020), 'scan 3' (updated and created on 04/08/2020), and 'scan 4' (updated and created on 04/08/2020). Each row has a circular 'View' button on the right. A yellow box highlights the entire table area. A large yellow box encloses the text 'Each Snapshot works like a separate Scan from this point on'. A yellow arrow points from this text to the 'View' button of the third snapshot ('scan 4'). Another yellow arrow points from the text 'Your new Snapshots are listed here' to the same 'View' button. A third yellow arrow points from the text 'Click on the eye to view' to the same 'View' button.


Each Snapshot works like a separate Scan from this point on

Your new Snapshots are listed here

Click on the eye to view

SCAN	UPDATED	CREATED
scan 2	04/08/2020	04/08/2020
scan 3	04/08/2020	04/08/2020
scan 4	04/08/2020	04/08/2020

5. Measurements

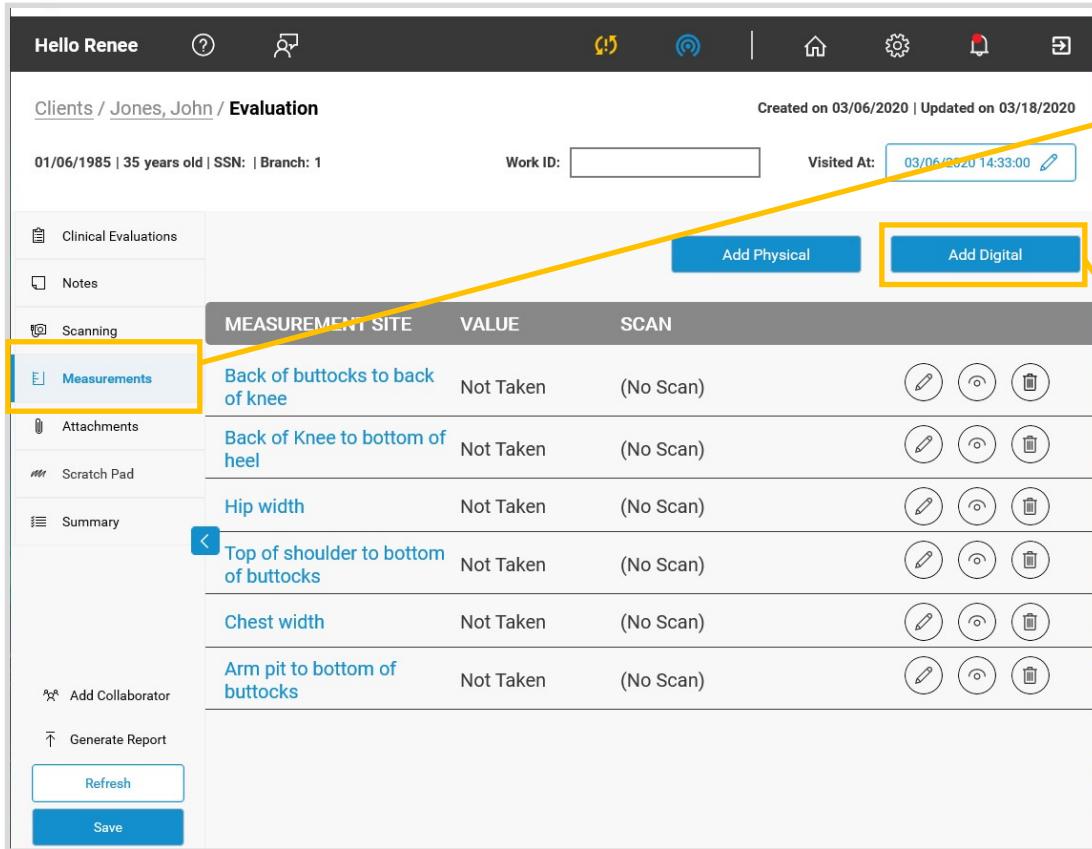


The screenshot shows the 'Clients / Jones, John / Evaluation' screen. On the left, a sidebar has 'Measurements' selected. At the top right, it says 'Created on 03/06/2020 | Updated on 03/18/2020'. Below that, there's a 'Work ID:' field and a 'Visited At:' field showing '03/06/2020 14:33:00'. Two blue buttons are at the top right: 'Add Physical' and 'Add Digital'. The main area is a table with columns: MEASUREMENT SITE, VALUE, and SCAN. It lists several measurements like 'Back of buttocks to back of knee', 'Back of Knee to bottom of heel', etc., all marked as 'Not Taken' with '(No Scan)' in the SCAN column. Each row has three icons: a pencil, a circular arrow, and a trash can.

Click on the Measurements Tab

Then click here to add a Digital Measurement

Add a Digital Point-To-Point Measurement

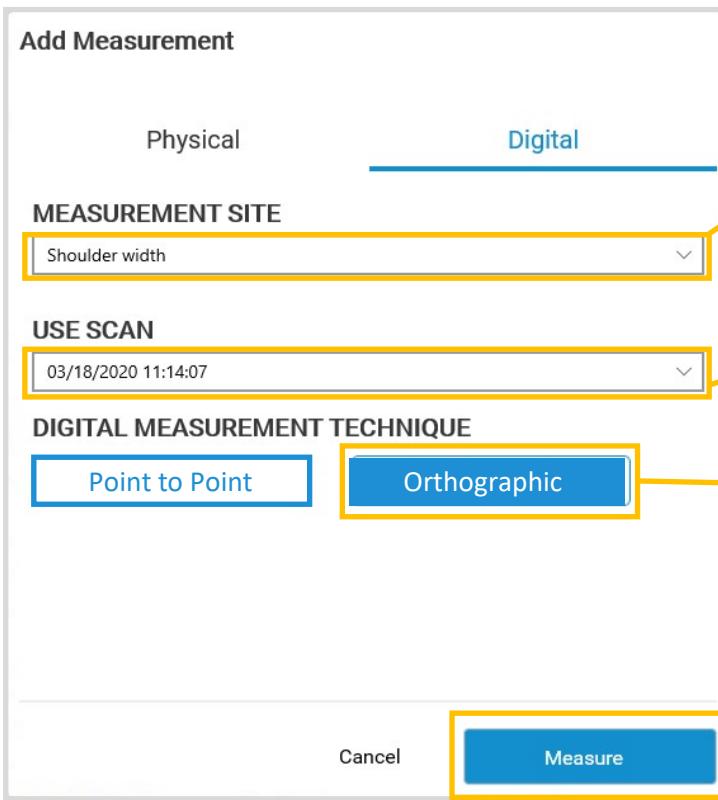

The dialog box is titled 'Add Measurement'. It has tabs for 'Physical' and 'Digital', with 'Digital' selected. Under 'MEASUREMENT SITE', a dropdown menu is open, showing 'Shoulder width'. Under 'USE SCAN', another dropdown menu is open, showing '03/18/2020 11:14:07'. Under 'DIGITAL MEASUREMENT TECHNIQUE', two buttons are shown: 'Point to Point' and 'Orthographic', with 'Orthographic' highlighted. At the bottom are 'Cancel' and a large blue 'Measure' button.


Select the measurement to be made

Select scan to be measured

Select Orthographic

Click on Measure


5. Measurements, continued . . .

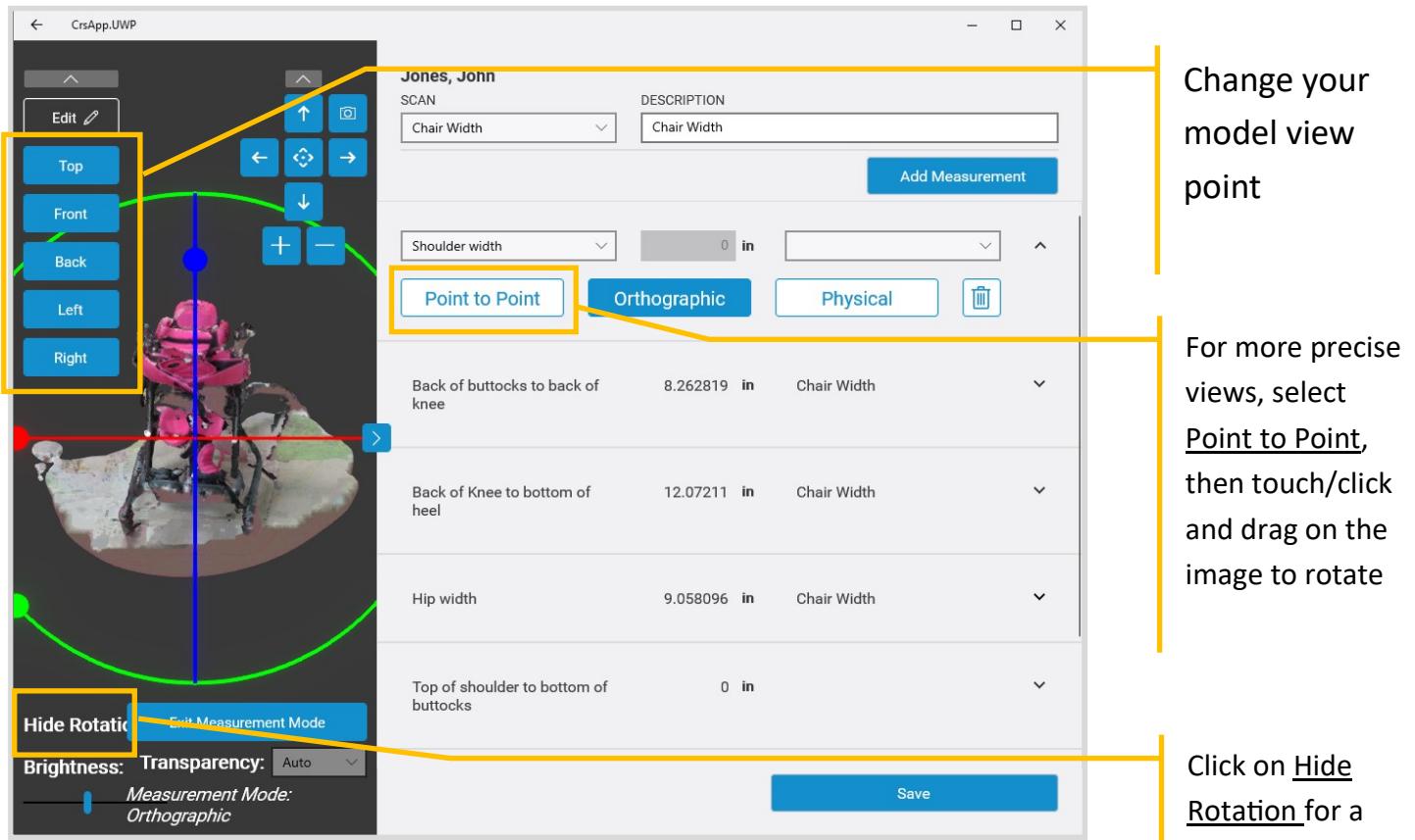


Click/touch to place measurement points on the image. Touch/click and drag to readjust

You will see the measurement displayed here

When you are happy with your measurement, click Save

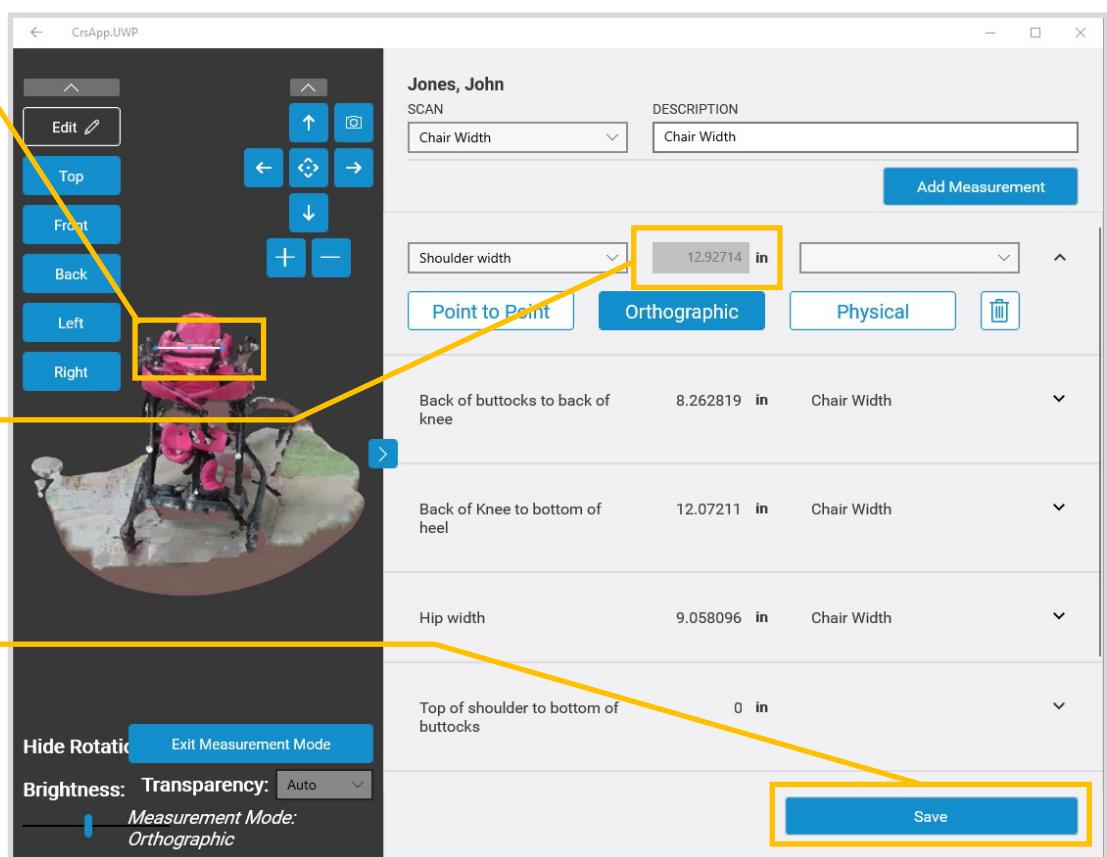


5. Measurements, continued . . .

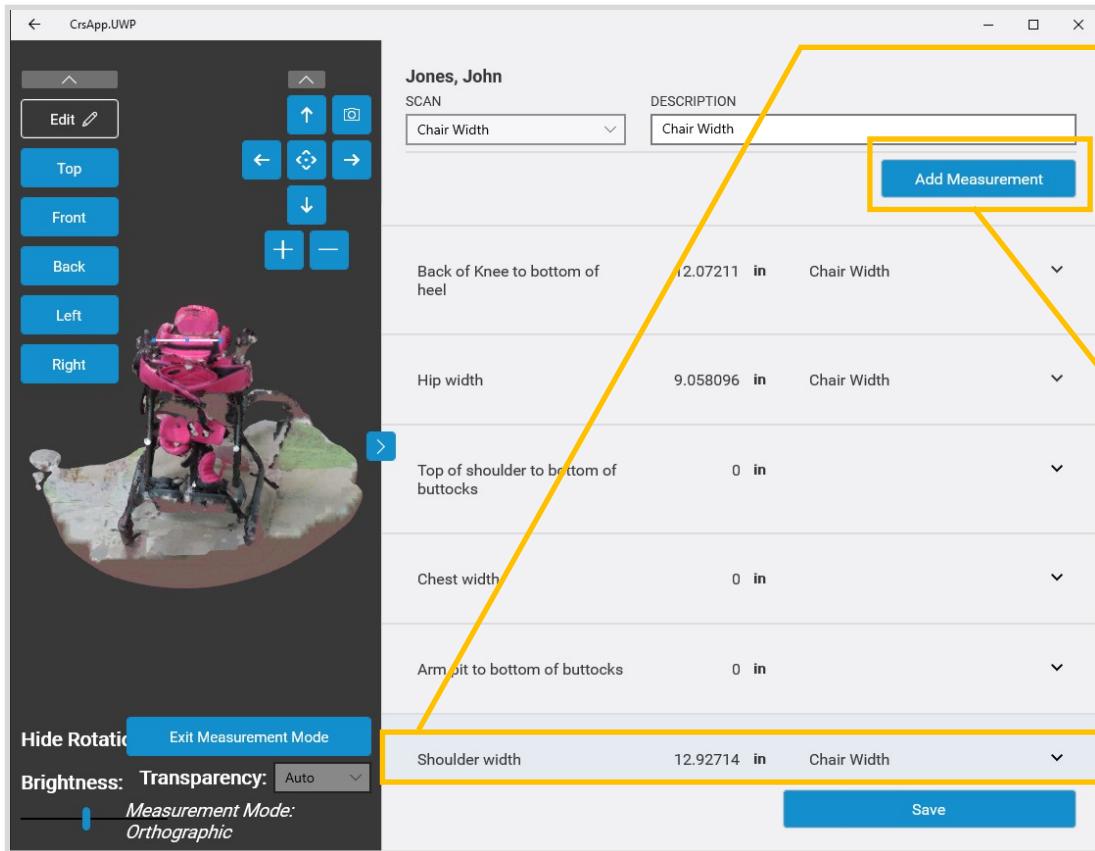

Your
measure-
ment will be
listed

Click to add
another
measurement

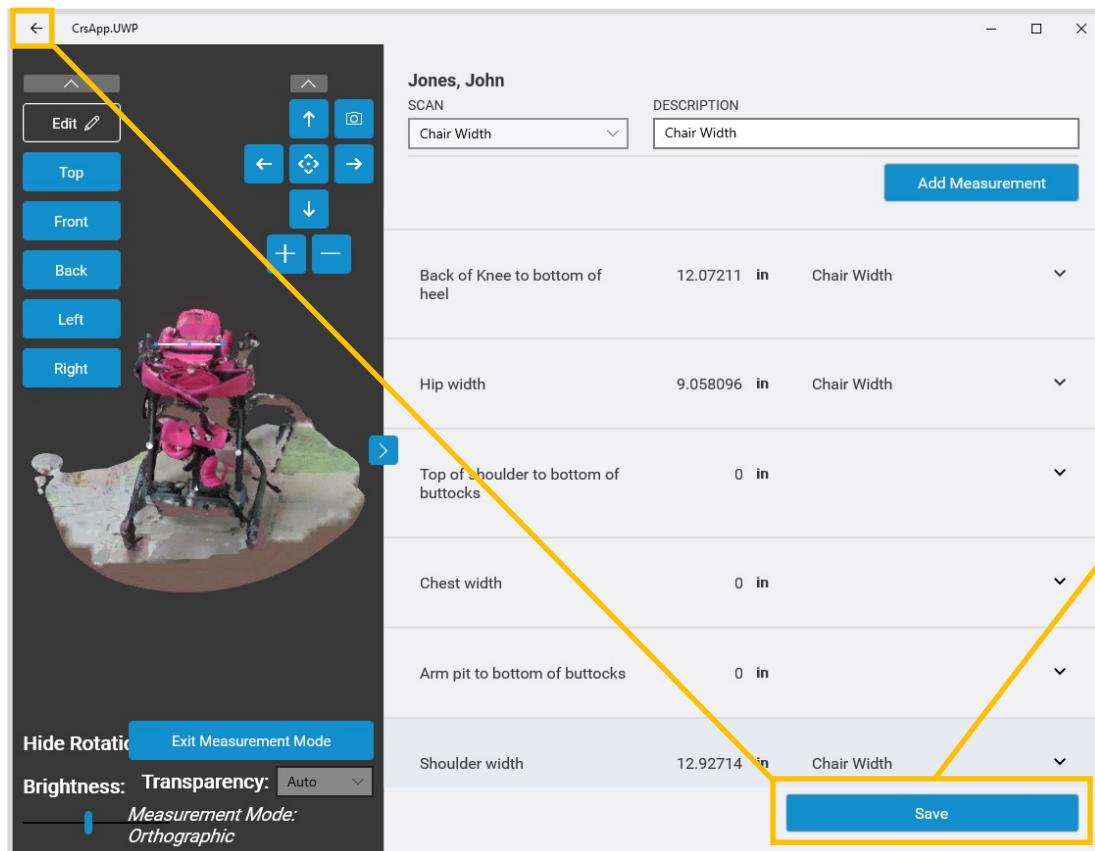

When you
are done taking
measurements,
click Save and
then the Back
Arrow

6. Summary Report

Clients / Jones, John / Evaluation

Created on 03/06/2020 | Updated on 03/20/2020

01/06/1985 | 35 years old | SSN: | Branch: 1

Work ID: _____ Visited At: 03/06/2020 14:33:00

MEASUREMENT SITE	VALUE	SCAN
Back of buttocks to back of knee	8.26 in	Chair Width <input type="button" value="Edit"/> <input type="button" value="Scan"/> <input type="button" value="Delete"/>
Back of Knee to bottom of heel	12.07 in	Chair Width <input type="button" value="Edit"/> <input type="button" value="Scan"/> <input type="button" value="Delete"/>
Hip width	9.06 in	Chair Width <input type="button" value="Edit"/> <input type="button" value="Scan"/> <input type="button" value="Delete"/>
Top of shoulder to bottom of buttocks	Not Taken	(No Scan) <input type="button" value="Edit"/> <input type="button" value="Scan"/> <input type="button" value="Delete"/>
Chest width	Not Taken	(No Scan) <input type="button" value="Edit"/> <input type="button" value="Scan"/> <input type="button" value="Delete"/>
Arm pit to bottom of buttocks	Not Taken	(No Scan) <input type="button" value="Edit"/> <input type="button" value="Scan"/> <input type="button" value="Delete"/>
Shoulder width	12.93 in	Chair Width <input type="button" value="Edit"/> <input type="button" value="Scan"/> <input type="button" value="Delete"/>
Top of head to bottom of buttocks	Not Taken	Chair Width <input type="button" value="Edit"/> <input type="button" value="Scan"/> <input type="button" value="Delete"/>

- Clinical Evaluations
- Notes
- Scanning
- Measurements
- Attachments
- Scratch Pad
- Summary

Click on
Summary

Clients / Jones, John / Evaluation

Created on 03/06/2020 | Updated on 03/20/2020

01/06/1985 | 35 years old | SSN: | Branch: 1

Work ID: _____ Visited At: 03/06/2020 14:33:00

Clinical Evaluations Notes Scanning Measurements Attachments Scratch Pad Summary

Add Collaborator

John Jones

01/06/1985 | 35 years old | SSN: | Branch: 1 Evaluation: 03/06/2020-14:33:00 UTC

Created On 03/06/2020 | Updated On 03/20/2020

Contacts

Name	GREG ORTIZ
Relation	Parent
Phone	8888675309
Email	greg@ortiz.com
Address	123 Cline Place, Montgomery AL

Measurements

Back of buttocks to back of knee 8.26 in

Toggle to include
Measurement
Images

Scroll down
to view

Click to
Generate a
Report



Technical Support

Need additional information or have questions?

There are several support resources available for users in the CRS app. Click the  icon in the header of the CRS application to access the following resources:

- Interactive Smart Search and Support Slides
- Video Tutorial Library
- Detailed Usage Guide
- Quick Start Guide

For additional help, our Tech Support Team can be reached during regular business hours by email or telephone:

Email CRS Technical Support at: crs.support@kimobility.com or call us at **(715) 997-9979**.

Our **technical support business hours** are:

Monday – Friday
8:00am – 5:00pm CST

Feedback

We want to hear from you!

If you have any feedback to share about the app, including ideas for new features, you can do so by clicking the  icon in the header of the CRS application.