

## CRS App Detailed Usage Guide

### Table of Contents

|                                       |           |
|---------------------------------------|-----------|
| <b>Overview</b>                       | <b>1</b>  |
| <b>CRS App Tour</b>                   | <b>1</b>  |
| <b>Detailed Usage</b>                 | <b>2</b>  |
| I: Logging In                         | 2         |
| II: Creating a New Client             | 3         |
| III: Creating an Evaluation           | 4         |
| IV: Capturing and Importing 3D Images | 6         |
| V: Orienting a 3D Scan                | 12        |
| VI: Recording Measurements            | 13        |
| VII: Attaching Images & Video Files   | 16        |
| VIII: Using Scratch Pad               | 17        |
| IX: Collaboration Features            | 18        |
| X: Summary View                       | 19        |
| <b>Technical Support</b>              | <b>20</b> |

---

## Overview

Complex Rehab Systems (CRS) provides software for Assistive Technology Professionals (ATPs) and Rehab Technology Suppliers (RTSs) who are working in the field of Complex Rehab. The CRS application allows users to collect detailed electronic records to assist in the design, documentation, review, and justification of custom mobility solutions. This includes the ability to capture a 3D scan of the client (the model) to "take them with you," and to act as a reference when designing the client's custom mobility solutions.

## CRS App Tour

### Home Screen

The home screen is where users can:

- Find existing client records
- Create new client records
- Add clients to the "my client" list by toggling the star shown to the right of the client's age.

### Client Screen

The client screen shows all data associated with a single client. This data includes:

- Client details
- A list of evaluations, where each evaluation contains detailed information about the client at a moment in time.

### Evaluation Screen

The evaluation screen houses all information related to a client's specific evaluation. This includes:

- Clinical evaluation form
- Notes
- 3D Scans
- Measurements
- Attached items (e.g. videos)
- Scratch pad – annotation of images or notes
- Summary


Users will also be able to select "Favorites" from the list of Evaluation cards.

Please Note:

The application automatically opens to the Client List (Home page) after logging in. Users starting from a different screen in any of the below steps can navigate to the Client list by clicking the Home page icon in the upper right menu panel.

## Detailed Usage

### I: Logging in

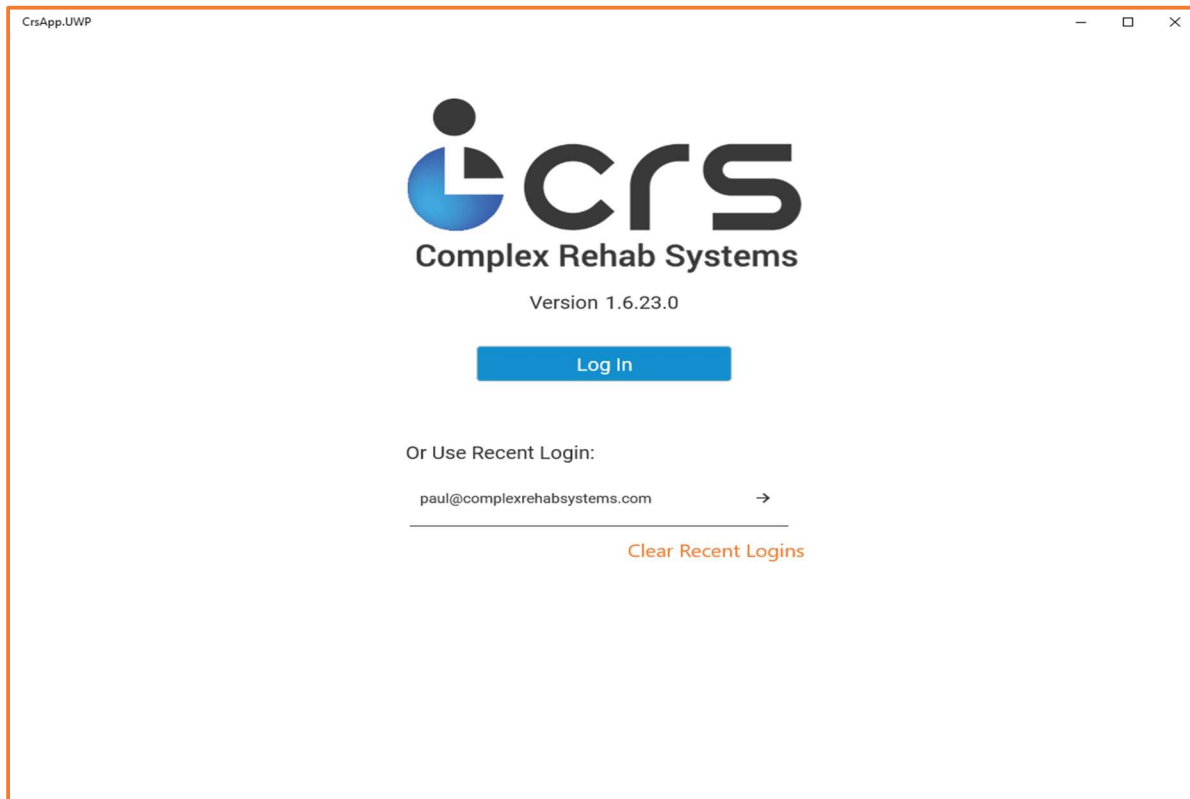

1. Launch the CRS application to get started.
2. Click the Log In button. **Note: If you have previously logged in, simply click the arrow next to your username.**
3. Enter your username.
4. Next, enter your password.
5. Click the Sign In button to proceed.

#### Offline login

- The first time a user logs in with an internet connection, they will be prompted to set up a PIN allowing for offline login.
- This PIN will be required at each login when no internet connection is detected.
- Users should configure a 4-6-digit PIN and click Confirm to proceed to the client dashboard.
- You may skip setting up a PIN, but you will be prompted each time you log in until one is created.

## II: Creating a New Client

CrsApp.UWP
— □ ×

Hello Renee
?
🗨
🔄
📶
|
🏠
⚙
🔔
📄

Clients

Search

🔍

Order by:

Last Updated A ▾

Branch:

Enter a Branch

Clients:

All clients ▾

🔄

Add Client

| NAME                     | UPDATED    | BRANCH | EXT. ID    | STATUS                 | SYNC       |
|--------------------------|------------|--------|------------|------------------------|------------|
| Hayes, Victor            | 04/10/2020 | 2      |            | Unsynchronized Changes | ☁          |
| NAME ↕                   | UPDATED ↕  | BRANCH | DOB        | AGE ↕                  | MY CLIENTS |
| Hayes, Victor            | 04/10/2020 | 2      | 02/03/1972 | 48                     | ☆ >        |
| Aardvark, Abner          | 04/10/2020 |        | 01/01/1970 | 50                     | ☆ >        |
| Zolo, Zane Z             | 04/10/2020 |        | 01/01/1970 | 50                     | ☆ >        |
| Demo, snapshot           | 04/10/2020 |        | 01/01/1970 | 50                     | ☆ >        |
| Jones, John              | 04/08/2020 | 1      | 01/06/1985 | 35                     | ★ >        |
| Jones, Jane              | 04/08/2020 | 1      | 01/06/1986 | 34                     | ★ >        |
| Kline, Jason Z           | 04/07/2020 | 2      | 02/01/1970 | 50                     | ☆ >        |
| Stevens, Renee Elizabeth | 03/20/2020 |        | 01/01/2000 | 20                     | ★ >        |


⏪
⏩

PAGE 1 OF 2

⏪
⏩

1. Click the 'Add Client' button.
2. Fill in client details and click Review. **Note: If a date of birth is not entered, the system will ask confirmation to leave the date of birth set at the default date.**
3. Review client details and click 'Confirm'. The application will look for possible duplicates and notify the user if any are found.
4. After confirmation users are taken directly to the Evaluation View.
5. If necessary, client details may be edited from the client screen.


### III: Creating an Evaluation

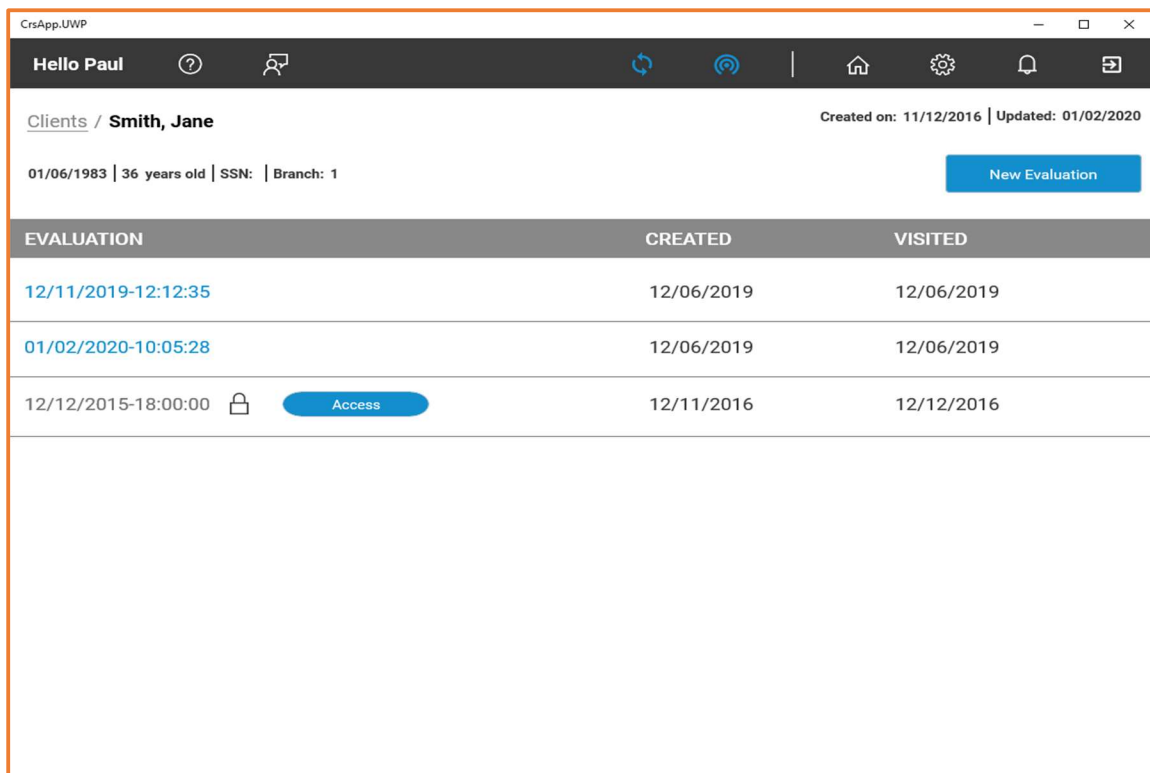

1. Navigate to a client view by selecting a client from the list.
2. Click the 'New Evaluation' button. This creates a blank evaluation record.
3. Open an existing evaluation by selecting the evaluation from the list in the client view.

Each evaluation has the following tabs:

- *Clinical Evaluations*: Enter contact information, home environments, medical history, and more.
- *Notes*: Add notes to the evaluation.
- *Scanning*: Capture new 3D scans or snapshots, view or import an existing 3D scan or snapshot.
- *Measurements*: View existing measurements and add new measurements to the evaluation.
- *Attachments*: Capture images and video or attach a file to the evaluation.
- *Scratch pad*: Annotate images or create and save additional notes.
- *Summary*: Summary view of all information entered in the evaluation.

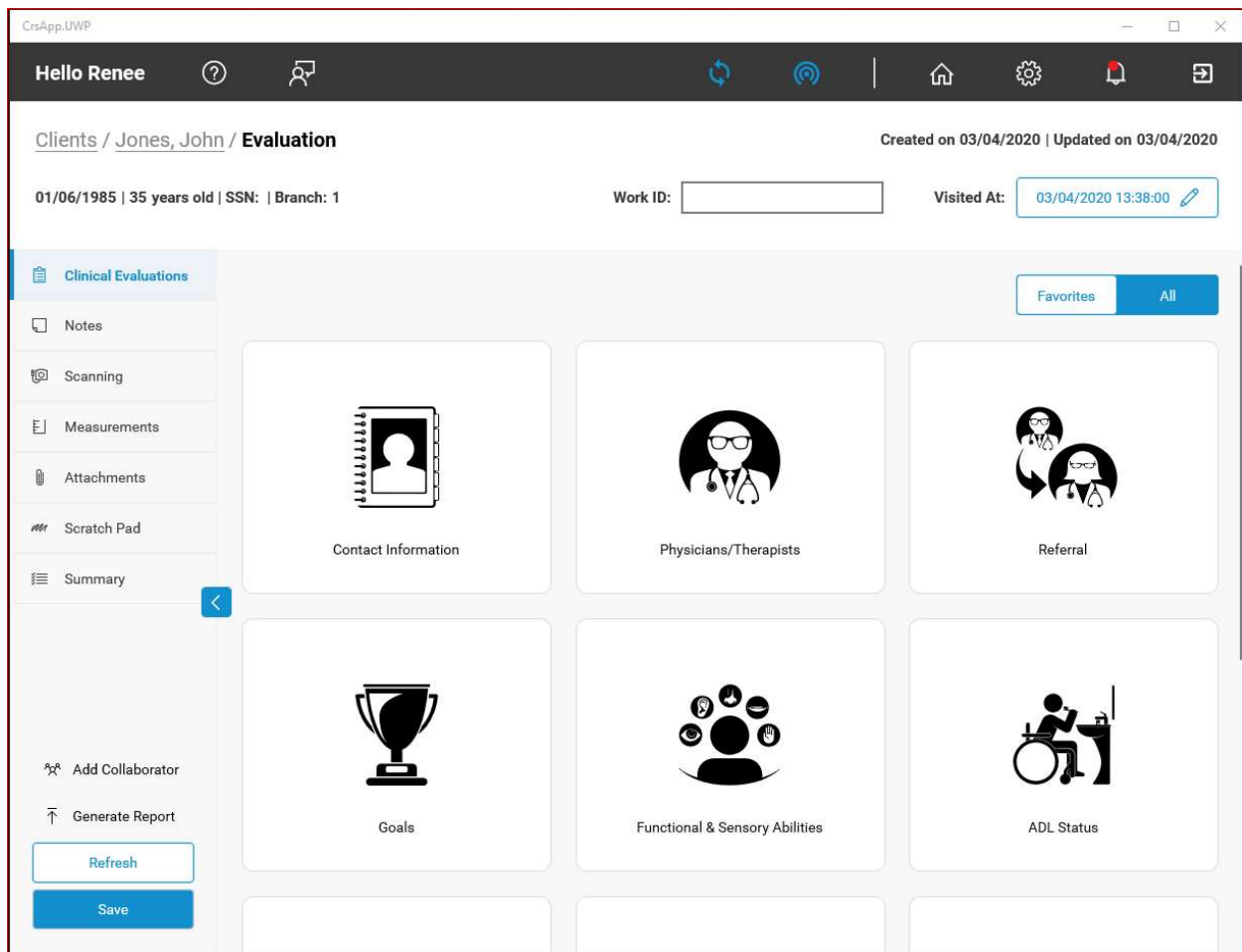

The following tabs are always available to the user during an evaluation:

- *Add Collaborator*: Allows a user to send a request to a colleague to work together on a specific evaluation.
- *Generate Report*: Creates a report with any information that has been entered during an evaluation.
- *Save*: Users are encouraged to save their work frequently and can do so at any time throughout the evaluation.

## IV: Capturing and Importing 3D Images

### Two Options: Scans and Snapshots

You have two choices when creating a 3D image of a client; Scan a Person or Object (Scans) and Take 3D Snapshots (Snapshots). Snapshots are like 3D pictures and are ideal for those who cannot remain still for 1-2 minutes. Scans creates a 3D full body image that can be viewed from any angle.

Go to page 10 for information on Snapshot.


Read on for more information on Scans.

### 3D Scans:

#### Best Practices for Capturing a Scan

To yield the best results, the following steps are recommended when taking a 3D scan in CRS App:

- Take the scan in a well-lit environment.
- The model must remain still during the entirety of the scan.
- Whenever possible, the model should wear light, colorful clothing as opposed to dark clothing. Tight-fitting clothes will also yield better results for accurate measurements.
- When setting the bounding box, all parts of the model must be completely inside the box.
- Maintain a consistent distance from the model throughout the scan.
- Limit the use of other background programs on your device during the scan.
- Ensure the sensor is aimed at the client and not blank space.
- Be careful to prevent abrupt movement of the sensor while scanning.

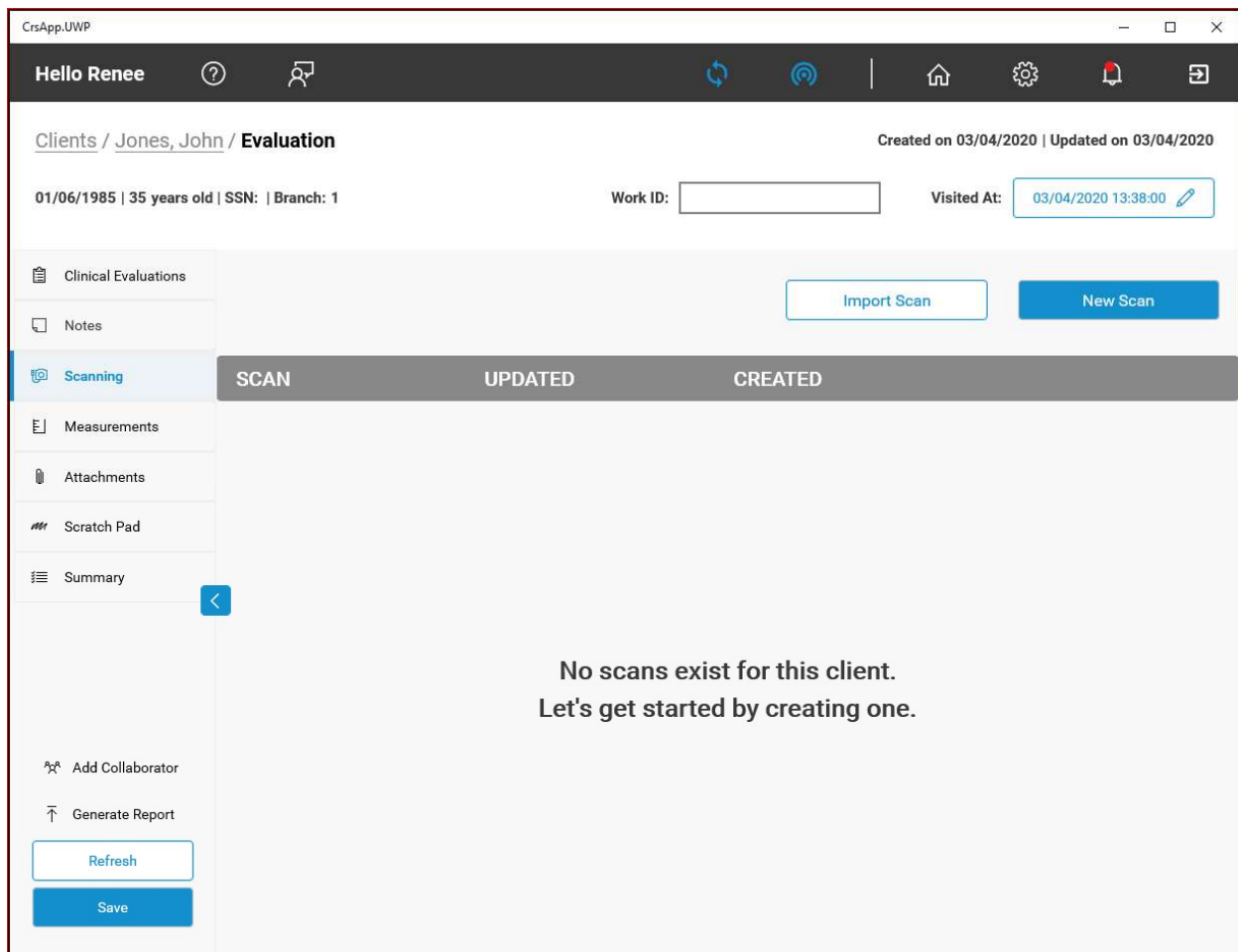


Multiple scans can be added to an evaluation record. These can be used for capturing linear measurements electronically (e.g. seat depth).

To capture a scan and attach it to the evaluation, the following steps are required:


- Navigate to a client evaluation and click the 'Scanning' tab.
- Launch the scanning application by clicking the 'New Scan' button.
- In the I want to... pop up, select Scan a Person or Object  
**Note: If a 3D sensor is not connected, the user will be prompted to connect a 3D Sensor.**
- Position the "bounding box" so it contains the entire model
  - Controls at the right allow the user to change the size and distance of the bounding box.
  - Note: Anything inside the bounding box will be included as part of the scan, anything outside the bounding box will be excluded from the scan.
  - Users should verify the bounding box contains the desired model by rotating the view.
  - The optimal distance between a user capturing a scan and the model is 2-4 feet. This distance should be maintained during the scan.



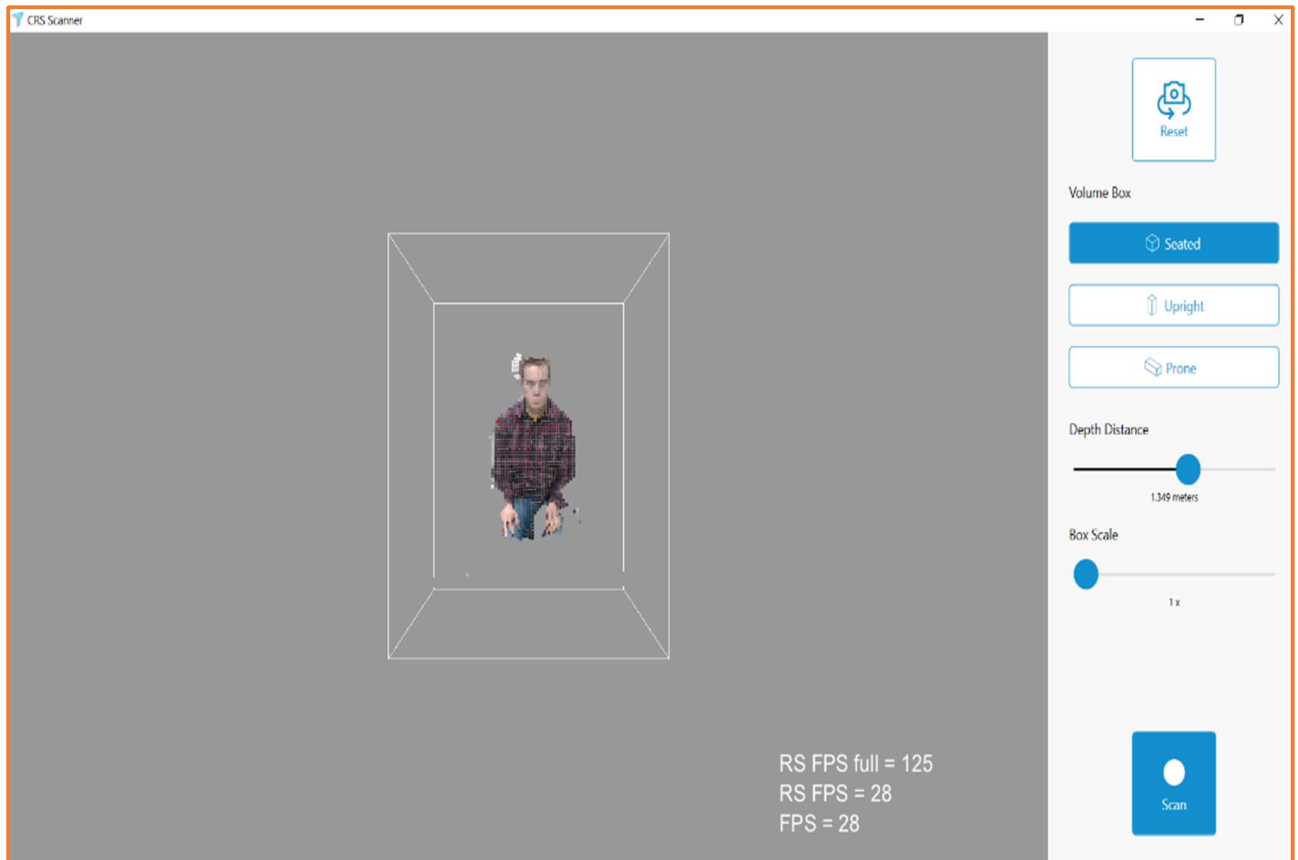


#### Bounding Box Controls:

- *Volume Box*: adjusts the shape of the bounding box.
- *Depth Distance*: moves the bounding box nearer to or farther from the sensor.
- *Box Scale*: increase/decrease the size of the bounding box.

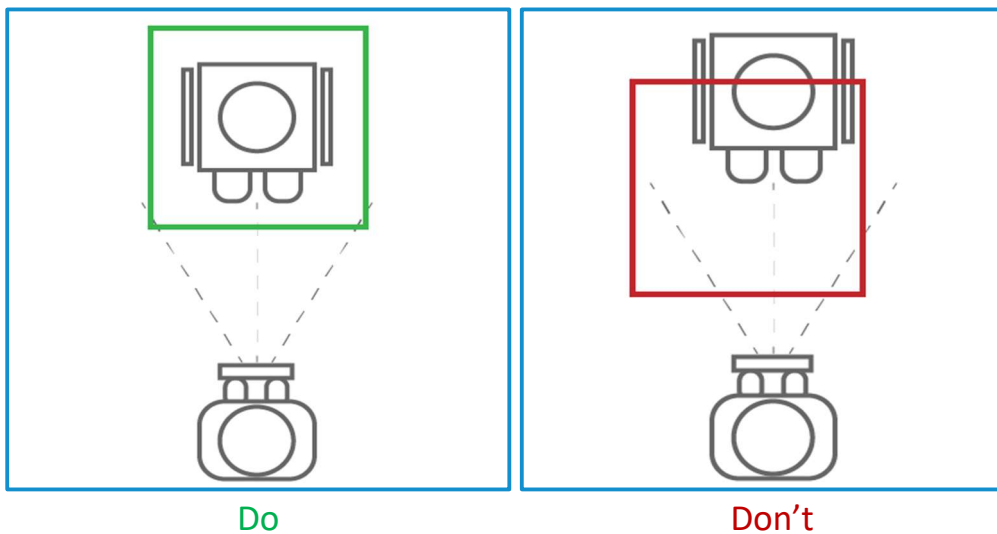

- Click the Scan button on the bottom right of the screen to start the scan
- Slowly walk around the client, being careful to move the sensor up and down to capture a complete model.

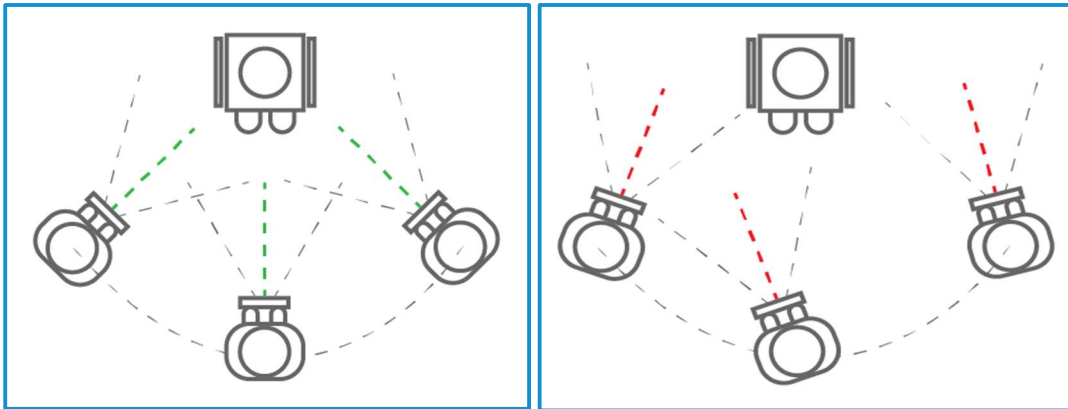

Do

Don't


- When satisfied with the model, click Complete. The data is now being formed into a 3D model and will appear on the device screen.
- If desired, the scan may be cropped prior to saving. When ready, click the 'Save' button to save the scan to the evaluation and return to the scanning tab.
- The 'eyeball' icon is used to open and view the scan, if desired.
- Save the evaluation now using the save button in the bottom left corner.
- A previously taken 3D scan can also be added to the evaluation by clicking the 'Import Scan' button.

### 3D Snapshots:

Snapshots are like a scan taken from only one perspective. Snapshots are taken instantaneously and are an ideal choice for the client that cannot remain still for 1-2 minutes. Taking a series of Snapshots from different directions and different areas of the body will give you a series of “scans” and each can be measured just like a regular scan.


To Take and Save a Series of 3D Snapshots:

1. While in an Evaluation, navigate to the Scanning tab and click on New Scan.
2. In the I want to... pop up, select Take 3D Snapshots
3. Point the scanning camera at the client and center them in your screen view.
4. Click on the Snapshot button to take a 3D Snapshot.
  - It is helpful to center and take Snapshots of areas of the body that include both points needed for a given measurement.

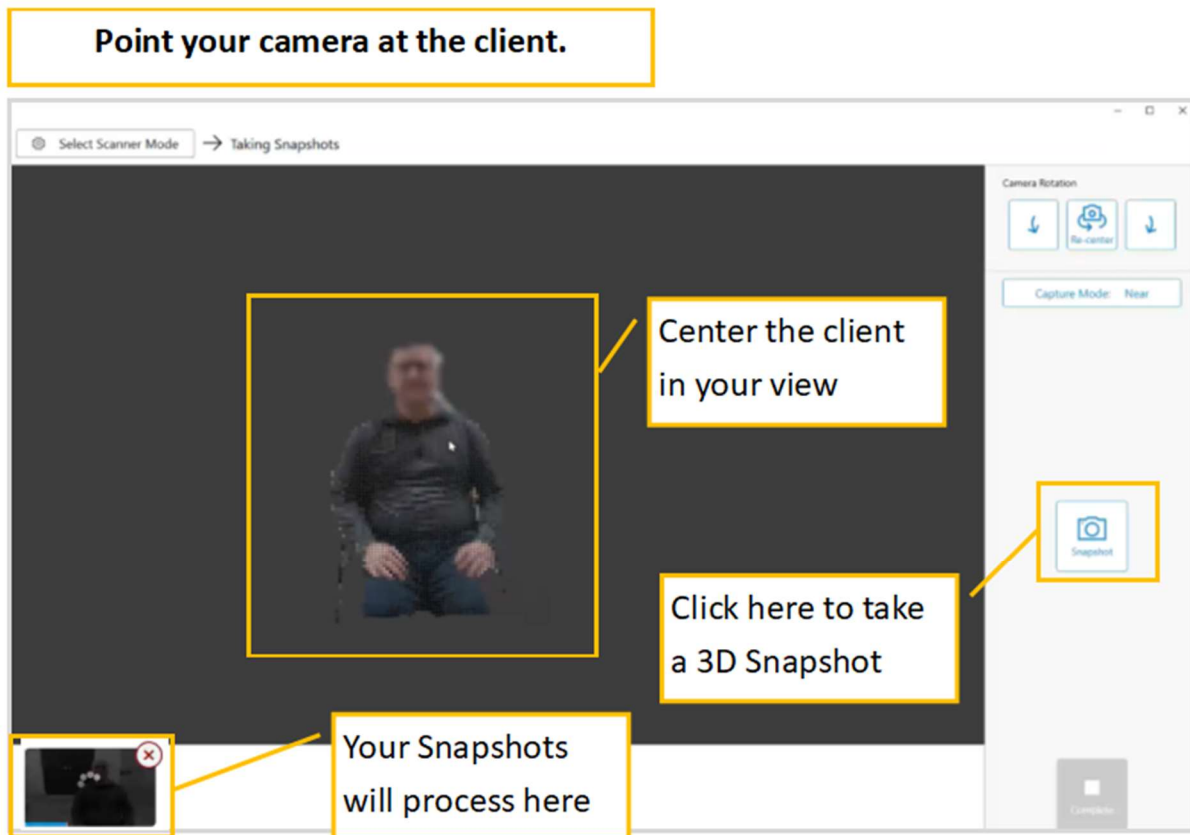

- The Snapshot will process in the lower left-hand corner of your screen.
5. Repeat steps 2-3 for each area needed for a measurement.
  6. When done, click on Complete.
    - It is recommended to click Complete after taking 6-8 Snapshots for faster processing times.

7. You will be taken to the Image Review screen. As they are done processing, click on each 3D image in the lower left-hand corner of your screen to view them.
  - You can rotate the 3D image left or right by clicking/touching next to the image and dragging left or right
8. If you are not satisfied with the 3D image, you can delete it by clicking on the red X in the thumbnail's right-hand corner. To delete all 3D images and start taking new Snapshots, click on Retake Scan and go back to step 2.

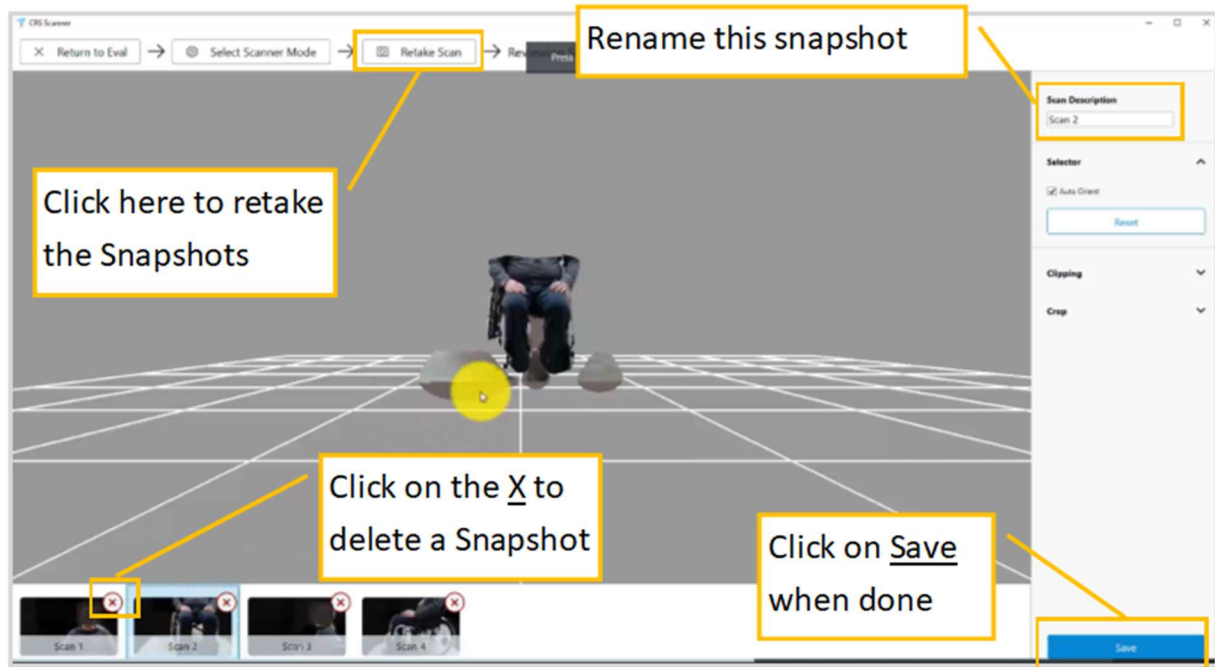


9. When you are satisfied with the 3D images, click on Save.
10. You will be returned to the Scanning screen. Each snapshot will be listed as a separate “scan”.
11. Click on the eye to the right of a listed scan to view it.
  - This is also a shortcut to taking Measurements.

Return to step 1 to take additional Snapshots, as needed.

## V: Orienting a 3D scan

All scans should auto-orient. If the scan does not auto-orient, or you wish to adjust the orientation of the scan, follow the instruction below.

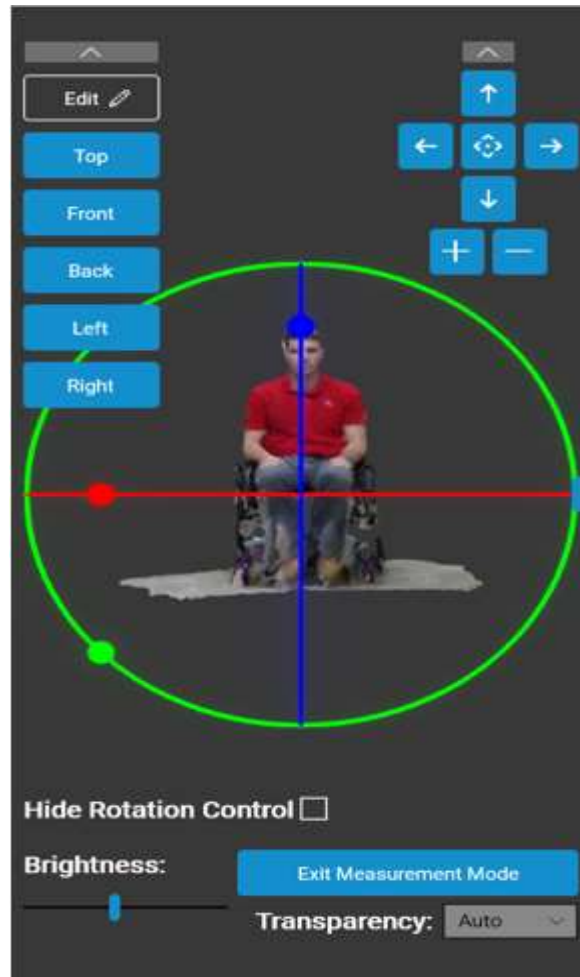

A scan can be aligned when it is open in the 3D scan viewer. To set the orientation of a scan:

1. Manually rotate the model on the screen such that it is viewed head on from the front. Click **'Edit'** (next to the Top/Front/Back/Left/Right view selection buttons) and click **'Set Front'**.  
When selecting any of these controls (e.g. Set Front) it will configure the model, so the current view is locked as that view (e.g. front). Only one of the views needs to be set; the remaining standard views will be calculated from there. CRS recommends always starting with the front view, and fine tuning the other views if necessary.
2. Once a view has been set, cycle through all the views (Front/Back/Left/Right/Top) by clicking on them individually and adjust alignment if necessary.
3. After the views are configured, users have the option to collapse the 'Set Views' controls by clicking the carrot above the controls.
4. Click the back arrow to return to the scanning tab.

## VI: Recording Measurements


The screenshot displays the 'CrsApp.UWP' interface. At the top, a header bar shows 'Hello Paul' and navigation icons. Below this, the breadcrumb 'Clients / Smith, Jane / Evaluation' is visible, along with the date 'Created on 12/06/2019 | Updated on 12/12/2019'. The client's details '01/05/1983 | 36 years old | SSN: | Branch: 1' are shown, followed by a 'Work ID' field and a 'Visited At' timestamp '12/06/2019 09:02:00'. A sidebar on the left contains navigation options: 'Clinical Evaluations', 'Notes', 'Scanning', 'Measurements' (selected), 'Attachments', 'Scratch Pad', and 'Summary'. The main area features a table with columns 'MEASUREMENT SITE', 'VALUE', and 'SCAN'. The table lists five measurements: 'Back of buttocks to back of knee', 'Arm pit to bottom of buttocks', 'Chest width', 'Top of shoulder to bottom of buttocks', and 'Hip width'. Each row shows 'Not Taken' for the value and '(No Scan)' for the scan, with edit, view, and delete icons. At the bottom, there are buttons for 'Add Physical', 'Add Digital', 'Add Collaborator', 'Generate Report', and 'Save'.

| MEASUREMENT SITE                      | VALUE     | SCAN      |
|---------------------------------------|-----------|-----------|
| Back of buttocks to back of knee      | Not Taken | (No Scan) |
| Arm pit to bottom of buttocks         | Not Taken | (No Scan) |
| Chest width                           | Not Taken | (No Scan) |
| Top of shoulder to bottom of buttocks | Not Taken | (No Scan) |
| Hip width                             | Not Taken | (No Scan) |

In the Measurements tab, users can add electronic records of measurements to the evaluation. There is no limit as to how many measurements can be added. Users have the option to create **digital** measurements (selecting *Point to Point* or *Orthographic*) or to enter **Physical** measurements (selecting *Physical*) that were previously taken using a measuring tape. Users may switch from physical measurement mode to digital by clicking the Physical or Digital tab within the Add New Measurement Screen.

To create a measurement from a scan, first determine the measurement type:

- **Point to Point** – this will measure on a straight line between any two points on the model. Users should place the start and end point of the measurement on the model after selecting this measurement mode. Point to point should be used when taking measurements of clients with asymmetry for accuracy; this measurement is bound directly to the model.
- **Orthographic** – the measurement will be projected on the model from one of the five standard views (Top/Front/Back/Right/Left). These measurements are not bound directly to the model and are instead snapped to the surface of a box surrounding the model.
- **Physical** – this measurement mode allows you to measure the physical person with measuring tape and enter the value into the text entry field.

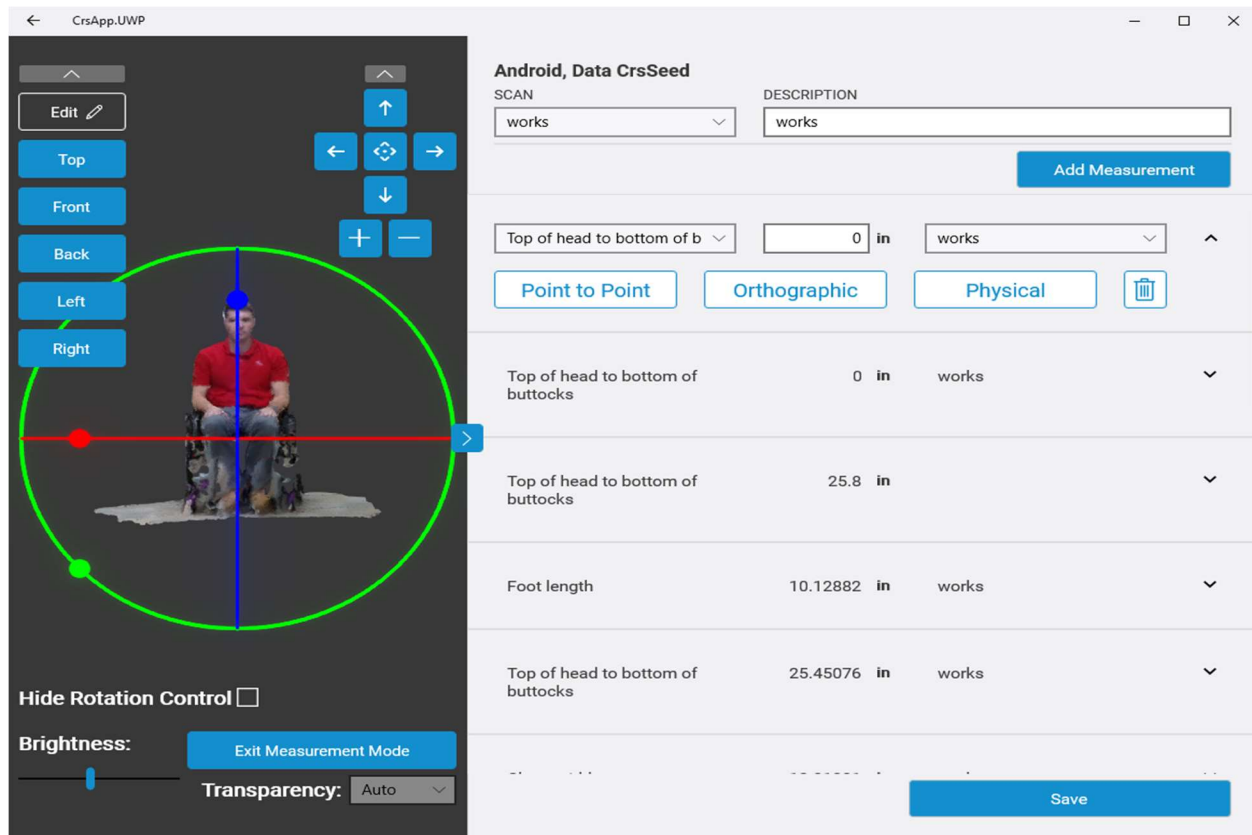

To Add new *Digital* measurements (Point to Point, Orthographic):

- Navigate to a client evaluation and click the Measurements tab.
- Click the 'Add Digital' button above the default measurement list.  
**Note: Each measurement must have a measurement site selected (e.g. "Shoulder Width") using the provided drop-down list.**
- Users must also select the desired 3D Scan that should be displayed for taking digital measurements.
- Select Point to Point or Orthographic mode (measurements cannot be taken if a mode has not been selected).
- Click Measure.

Both point-to-point and orthographic measurement types are taken by clicking two points on the viewer. Once both points have been placed, the measurement will populate the value of the corresponding measurement in the measurement list. These are both digital measurements and can be edited modifying the blue control points in the 3D viewer.

To Add new *Physical* measurements:

- Navigate to a client evaluation and click the Measurements tab.
- Click the 'Add Physical' button above the default measurement list.

**Note: Each measurement must have a measurement site selected (e.g. “Shoulder Width”) using the provided drop-down list.**


- Enter the measurement value and click Save. Each measurement must have a value before it can be saved.

*Note:*

- Measurements that have been entered manually can be edited using the text entry field.
- Measurements that have been taken digitally may not be edited directly using the text field; these must be edited using the control points in the 3D scan viewer.
- Measurements can be removed by clicking the ‘trashcan’ icon at the right side of the measurement within scan view or on the measurement tab.
- If the measurement control points fall inside the model, the viewer will become transparent to ensure users can still edit the measurements properly.
- After clicking ‘**Save**’ to save the measurements, use the back arrow in the upper left corner to go back to the measurements tab.



## VII: Attaching Images and Video

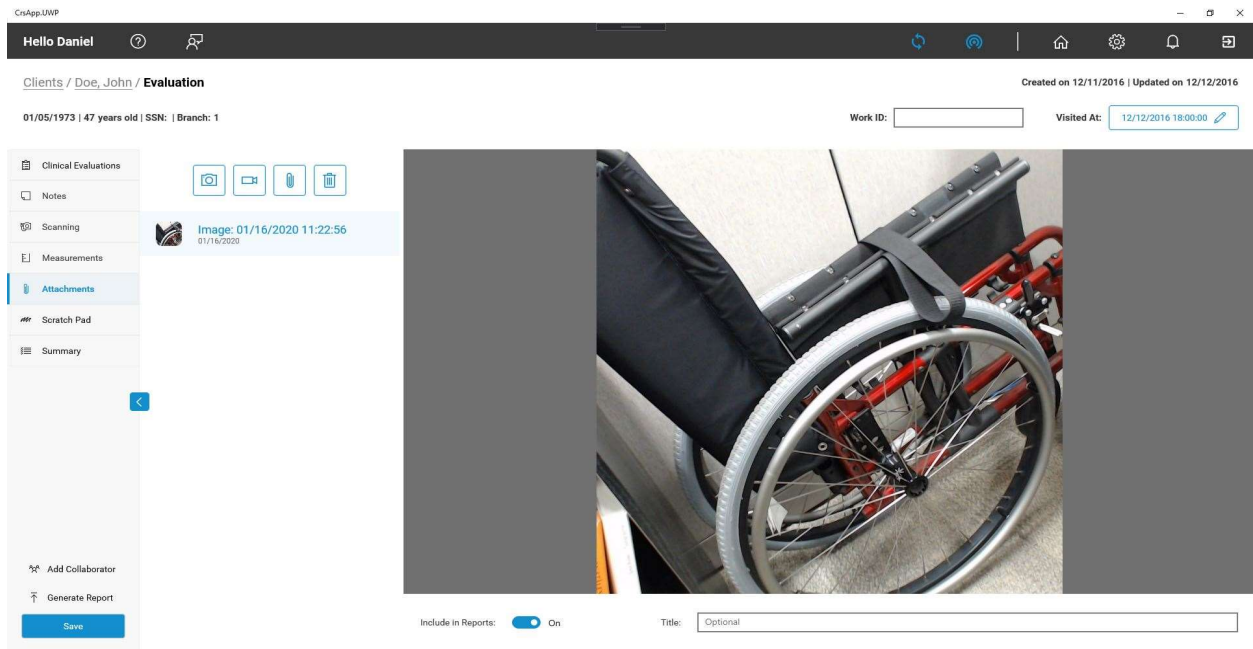

Videos and images can be attached directly to an evaluation on the Attachments tab. They are captured using standard device camera utilities and are automatically attached to the evaluation.


1. Navigate to a client evaluation and click the Attachments tab.
2. Select the desired media to attach (e.g. photo, video, etc.).
3. When choosing photo or video, cycle through the camera utilities (front camera, rear camera, etc.) using the built-in windows camera application.
4. For photo, click the camera icon to take a picture and then select 'Done' to save the new image to the attachments tab.
5. For videos, use the start/end record buttons and then select 'Ok' to save.

### Note:

- Images/video that have already been captured can be viewed by selecting them from the list provided.
- Images and videos have a preview screen; attaching other file types (e.g. PDF) will open in the device default viewer.
- Users have the option to type a description/title regarding the attached item in the bottom right corner after selecting it from the list.
- Users may also choose to exclude attachments in PDF reports by switching the toggle (next to the description/title box) to the 'Off' position. The application defaults to the 'On' position.
- Click 'Save' to update.

## VIII: Using Scratch Pad

1. Navigate to a client evaluation and click the 'Scratch Pad' tab.
2. Click 'New Sheet' to create a new scratch pad sheet.
3. Entering a title in the upper left corner is optional.

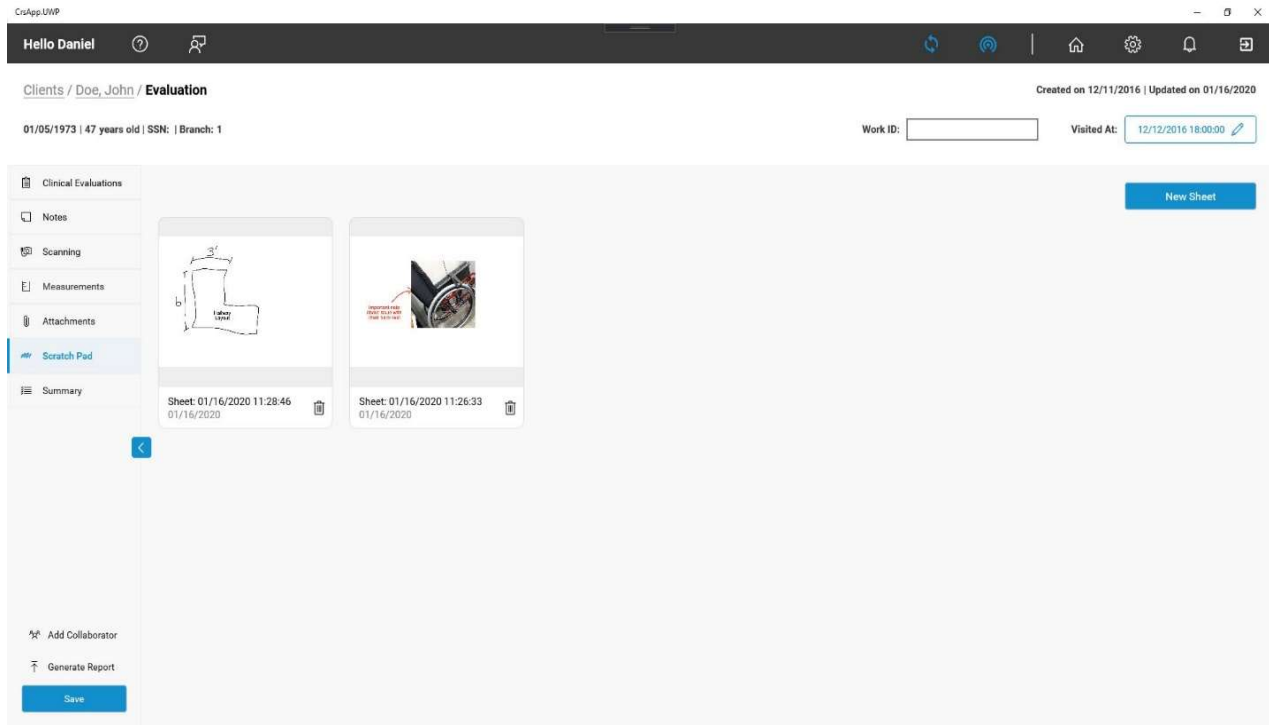


Each sheet will include the following drawing tools:

- *Selector*: Select any of the existing drawing elements that have already been added to the sheet
  - *Pen*: Freehand pen for easy annotation
  - *Text Entry*: Add text to the sheet
  - *Picture*: Add images to the sheet
4. Use the drawing tools to attach/annotate images on the sheet.
  5. Save the sheet to the evaluation by clicking the 'Save' button in the bottom right corner.
  6. Click the back arrow in the upper left corner to go back to the scratch pad tab after saving.

## IX: Collaboration Features

There are multiple ways for users to collaborate with each other on evaluations. Want to give another user access to an evaluation?

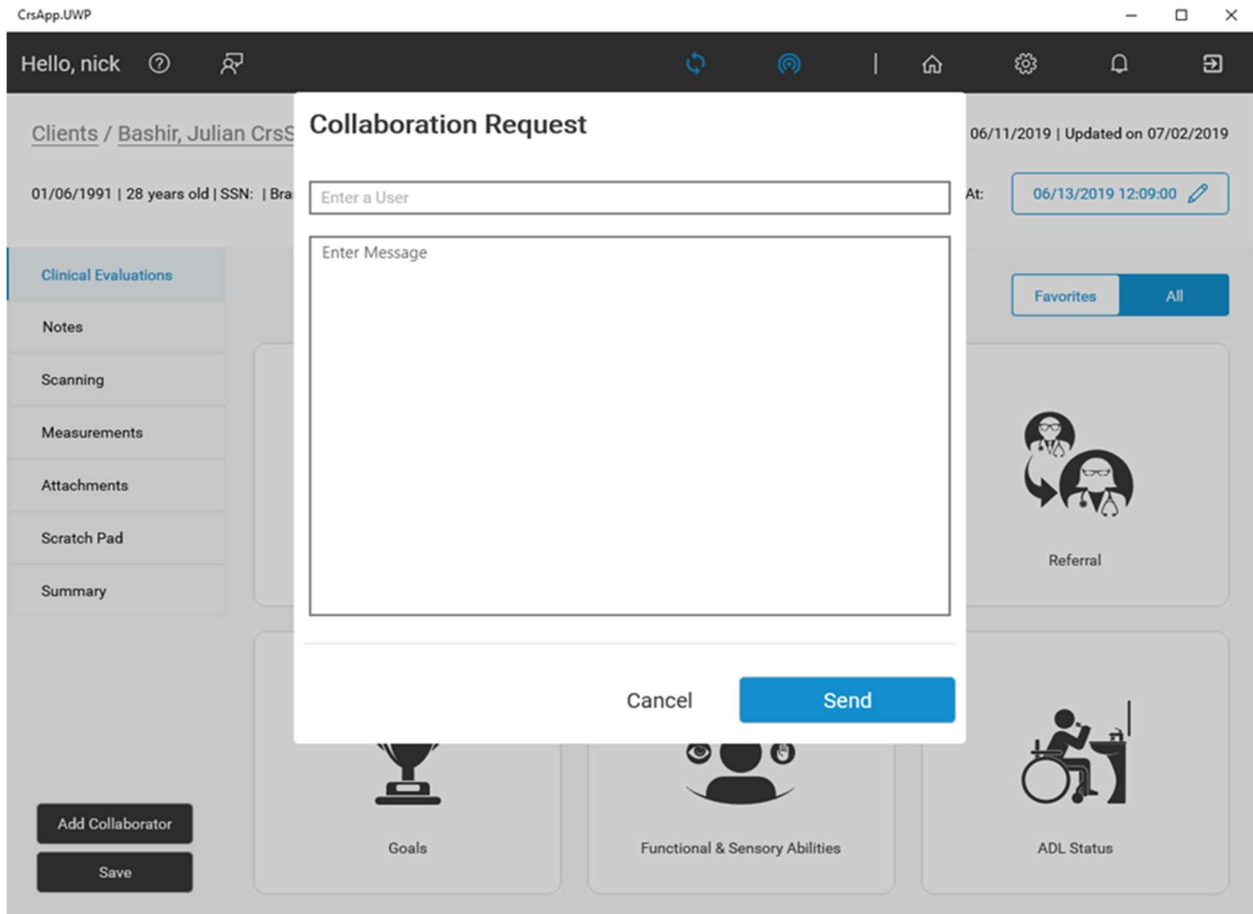

Complete these steps to collaborate on an evaluation:

- Click the 'Add Collaborator' button in the bottom left of the screen when working on an evaluation. This will bring up the collaboration request window.
- Enter the name of your colleague you would like to collaborate with on an evaluation.
- Add an optional message in the dialogue box and click 'Send'.
- The recipient of the collaboration request will receive an email with instructions on how to access the specific evaluation.

See an evaluation you are locked out of? A locked evaluation is demonstrated by a padlock icon and an 'Access' button. Complete these steps to request access to a locked evaluation:


- Click the 'Access' button.
- A new window will appear with space for a custom optional message to be sent with the access request.
- Click the 'Request Access' button to send the request.


- The lead ATP on the evaluation will receive an email asking them to approve or deny the access request.
- Once a user has access to an evaluation, they will have full read/write access.

## X: Summary View

At any point in the evaluation users may click 'Summary' to:

- View a summary of the client evaluation,
- Export a report,
- or make any additional notes.





01/06/1973 | 46 years old | SSN: | Branch: 1  
Created On 12/12/2016 | Updated On 12/06/2019

Evaluation: 12/12/2016-18:00:00 UTC  
Work Id:

### Contacts

|          |                       |
|----------|-----------------------|
| Name:    | MARTY MCFLY           |
| Relation | Self                  |
| Email    | mmcfly@backintime.com |

### Measurements

|                                   |          |
|-----------------------------------|----------|
| Back of buttocks to back of knee  | 23.36 in |
| Top of head to bottom of buttocks | 32.22 in |
| left shoulder to neck             | 5.7 in   |
| Head width                        | 6.02 in  |


### Residences

|               |                         |
|---------------|-------------------------|
| Location Type | TERTIARY RESIDENCE      |
| Home Type     | LONG-TERM CARE FACILITY |



## Technical Support

### Need additional information or have questions?

There are several support resources available for users in the CRS app. Click the  icon in the header of the CRS application to access the following resources:

- Interactive Smart Search and Support Slides
- Video Tutorial Library
- Detailed Usage Guide
- Quick Start Guide

For additional help, our Tech Support Team can be reached during regular business hours by email or telephone:

Email CRS Technical Support at: [crs.support@kimobility.com](mailto:crs.support@kimobility.com) or call us at **(715) 997-9979**.

Our **technical support business hours** are:

Monday – Friday

8:00am – 5:00pm CST

## Feedback

### We want to hear from you!

If you have any feedback to share about the app, including ideas for new features, you can do so by clicking

the  icon in the header of the CRS application.